


MPS4264

MINIATURE PRESSURE SCANNER

HARDWARE AND SOFTWARE MANUAL

SOFTWARE VERSION 2.07

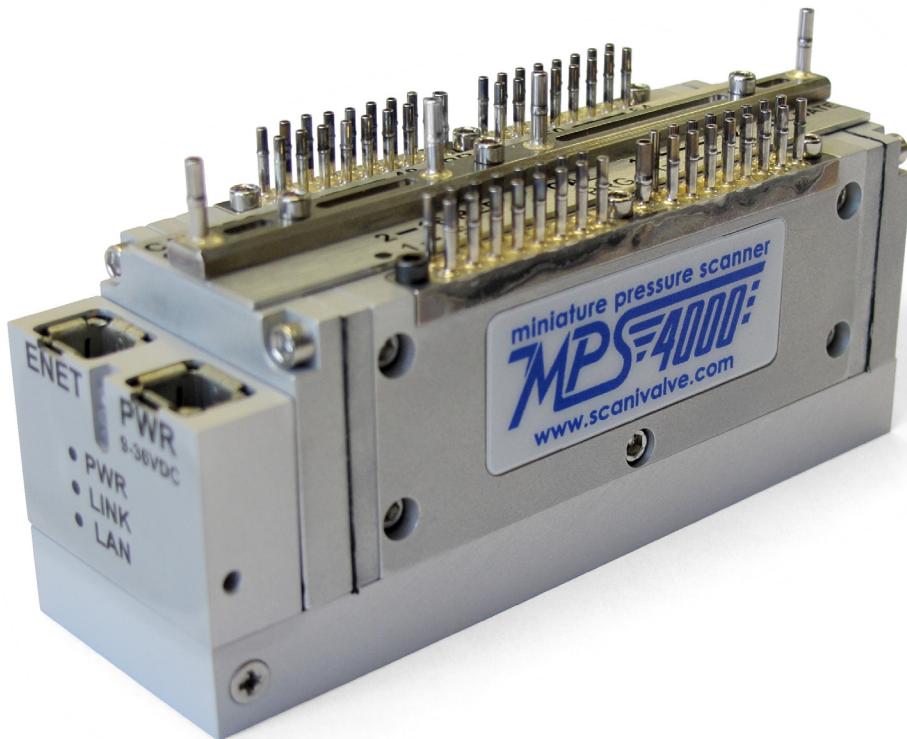

Scanivalve

PREFACE

WARNINGS, CAUTIONS AND NOTES



The WARNING! symbol indicates that danger of injury for persons and the environment and/or considerable damage (mortal danger, danger of injury) will occur if the respective safety precautions are not taken.



The CAUTION ! symbol indicates danger for the system and material if the respective safety precautions are not taken.



The ESD note symbol indicates that proper precautions for handling Electrostatic Sensitive Devices needs to be taken when performing the related operation. This includes the use of grounded work surfaces and personal wrist straps to prevent damage to sensitive electronic components.

WARRANTY

Scanivalve Corporation, Liberty Lake, Washington, hereafter referred to as Seller, warrants to the Buyer and the first end user that its products will be free from defects in workmanship and material for a period of twelve (12) months from date of delivery. Written notice of any claimed defect must be received by Seller within thirty (30) days after such defect is first discovered. The claimed defective product must be returned by prepaid transportation to Seller within ninety (90) days after the defect is first discovered. Seller's obligations under this Warranty are limited to repairing or replacing, at its option, any product or component part thereof that is proven to be other than as herein warranted.

Surface transportation charges covering any repaired or replacement product or component part shall be at Seller's expense; however, inspection, testing and return transportation charges covering any product or component part returned and redelivered, which proves

not to be defective, shall be at the expense of Buyer or the end user, whichever has returned such product or component part.

This Warranty does not extend to any Seller product or component part thereof which has been subjected to misuse, accident or improper installation, maintenance or application; or to any product or component part thereof which has been repaired or altered outside of Seller's facilities unless authorized in writing by Seller, or unless such installation, repair or alteration is performed by Seller; or to any labor charges whatsoever, whether for removal and/or reinstallation of the defective product or component part or otherwise, except for Seller's labor charges for repair or replacement in accordance with the Warranty. Any repaired or replacement product or component part thereof provided by Seller under this Warranty shall, upon redelivery to Buyer, be warranted for the unexpired portion of the original product warranty.

THIS WARRANTY IS IN LIEU OF AND EXCLUDES ALL OTHER WARRANTIES, EXPRESS OR IMPLIED, ARISING BY

OPERATION OF LAW OR OTHERWISE, INCLUDING THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, AND IN NO EVENT SHALL SELLER BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES.

In the event of a failure:

- 1) Notify Scanivalve Corporation, Customer Service Department. Include model number and serial number. On receipt of this information, service data or shipping instructions will be forwarded. This may be transacted by telephone or e-mail.
- 2) On receipt of shipping instructions, forward the product, transportation prepaid. Repairs will be made and the product returned.
- 3) All shipments should be made via "Best Way". The product should be shipped in the original packing container or wrapped in protective material and surrounded by a minimum of four (4) inches of a shock absorbing material.

TRADEMARKS ® AND COPYRIGHTS ©

Scanivalve is a registered trademark of Scanivalve Corporation.

All other brand and product names are trademarks or registered trademarks of their respective companies.

PACKAGING FOR SHIPMENT

If the product must be shipped, whether being returned to Scanivalve or relocated to another location it must be packaged properly to minimize the risk of damage. The recommended method of packing is to place the instrument in a container, surrounded on all sides with at least four inches of shock attenuating material such as Styrofoam peanuts.

IMPORTANT NOTICE

Please note that the product specifications and other information contained in this manual are subject to change without notice. Scanivalve Corporation makes an effort and strives to provide complete and current information for the proper use of the equipment. If there are any questions regarding this manual or the proper use of the equipment, contact Scanivalve Corporation.

CONTACT INFORMATION

If there are any questions or concerns regarding any Scanivalve product please do not hesitate to contact us at the following:

Scanivalve Corp.
1722 N. Madson Street
Liberty Lake, WA 99019
Telephone: (800)935-5151 (509)891-9970
Fax: (509)891-9481
scanco@scanivalve.com
www.scanivalve.com

Scanivalve Corporation is an ISO 9001:2015 certified company.

[INTENTIONALLY LEFT BLANK]

TABLE OF CONTENTS

PREFACE	IV
Warnings, Cautions and Notes	iv
Warranty	iv
Trademarks ® and Copyrights ©	v
Packaging for Shipment	v
Important Notice	v
Contact Information	v
SECTION 1: SPECIFICATIONS	7
General Specifications	7
Environment Specifications	7
SECTION 2: INTRODUCTION	9
Overview	9
Module Configurations	9
Pressure Range	9
Valve Actuation	9
Pressure Sensors	9
Pneumatic Valve	10
Pneumatic Inputs	11
Power & Ethernet Connections	11
SECTION 3: OPERATION	13
Unpack & Inventory	13
Mounting	13
Environmental Considerations	13
Warm-up	14
Communications	14
Serial Communications	14
Ethernet Communications	14
PC - TCP/IP	15
PC - UDP	15
PC - ScanTel	15
PC - LabVIEW	15
PC - HyperTerminal	15
FTP	15
Setting Up an Ethernet Connection	16
Using a Static IP address	16
Using a Dynamic IP address	16
Web Server: Connecting	17
Web Server: Configuring Scan Settings & Scanning Data to a File	17
Web Server: Scanning Data to the Screen	18
Web Server: Calibrations & Validations	19
Zero Cal	20

Web Server: "Remote" Screen	20
Web Server: "Terminal" Screen	21
FTP Scan Data Transfer	21
The Binary TCP Server	21
Commands	21
Scan Data	22
LabVIEW Scan Data	22
LabVIEW Data Buffering	22
Example LabVIEW	22
Multiple Unit Operations	22
Multiple Unit Scan Command	23
Multiple Unit Stop Command	23
Multiple Unit Find Command	23
Scanning with an External Trigger	23
Setting up an External Frame Trigger	23
Setting up an External Scan Trigger	23
Auto Scan Option	24
Statistical Scan Option	24
Fast Scan Option (850Hz - 2500Hz)	24

SECTION 4: HARDWARE**27**

Pneumatic Connections	27
Basics of Connecting Input Tubes	27
Pressure Measurement (Px) Ports	28
Calibration (CAL) Port	28
Reference (REF) Port	28
Purge (PRG) Ports	28
Calibrate Mode Control (CAL-PRG CTL) Port	28
Measurement Mode Control (Px CTL) Port	28
Removable Pressure Input Headers	29
Calibration Valve	29
Electrical Connectors	30
Power Connector & Power Requirements	30
Ethernet Connector	30
MPSTCU-Thermal Control Unit	30

SECTION 5: SOFTWARE**33**

Overview	33
Command Format	33
Output Scan Data Formatting and Routing	33
ASCII Data Formats	35
Basic ASCII Data Output	35
Formatted ASCII Data Output	36
ASCII CSV Data Output	36
Binary Data Formats	37
Binary Data Format	37
Statistical Binary Data Format	38
Fast Scan Binary Data Format	39
File Names and Formats	40

General Control Commands	41
Boot Loader Version	41
Firmware Files Check	41
Format Flash Disk	42
Multiple Unit Find	42
Multiple Unit Stop	42
Reboot	42
Status	43
Stop	43
TAR File Unpack	43
Valve State	44
Version	44
PTP (Precision Time Protocol) Commands	45
Clear Accumulated PTP Statistical Data	45
Display PTP Statistics	45
Get PTP Time	45
Get Scan Start Time	45
Get UTC Offset	46
Set PTP Time	46
Scan and Calibration Commands	47
Calibration	47
Multiple Unit Calibration	47
Calibration & Validation	47
Multiple Unit Calibration & Validation	48
Max-Min Calibration Table Values	49
Read Current Temperature	49
Scan	50
Multiple Unit Scan	50
Software Frame Trigger	50
Validation	51
Multiple Unit Validation	51
Zero Offset Calibration	51
Save Commands	52
Save	52
Scan Variables	53
List Scan Settings (LIST S)	53
Set Rate (RATE)	53
Set Frames Per Scan (FPS)	53
Set Scan Units (UNITS)	54
Set Scanning Data Format (FORMAT)	55
Set Trigger (TRIG)	55
Set Enable FTP (ENFTP)	55
Set Options (OPTIONS)	56
UDP Variables	57
List UDP Settings (LIST UDP)	57
Set UDP Enable (ENUDP)	57
Set the UDP Target Address and Port (IPUDP)	57
Identification Variables	58
List Identification Settings (LIST ID)	58
Set Serial Number (SN)	58
Set Nominal Pressure Range (NPR)	58
Set the Multiple Device Multicast Address (MCAST)	58

IP Variables	59
List IP Settings (LIST IP)	59
Set IP Address of the MPS (IPADD)	59
Set Subnet of the MPS (SUBNET)	59
Set Login Name (LOGIN)	60
Set Password for Login Name (PASSWORD)	60
Set Login Name 1 (LOGIN1)	60
Set Password for Login Name 1 (PASSWORD1)	60
Set Allow Anonymous FTP Login (ALLOWANON)	60
Set Application File Name (APP)	61
Set Gateway (GW)	61
Set DNS Name (NAME)	61
Precision Time Protocol Variables (LIST PTP)	62
LIST Precision Time Protocol Settings (LIST PTP)	62
Set PTP Enable (PTPEN)	62
Set PTP Statistical Output (STAT)	62
Set Scan Start Time (SST)	63
Set Scan Start Date (SSD)	63
Set UTC Offset (UTCOFFSET)	63
Conversion Table Variables	64
List Conversion Table Settings (LIST T)	64
Set Conversion Table Point (PT)	64
Calibration Variables (LIST C)	65
List Calibration Settings (LIST C)	65
Set Number of Points in Conversion Table (NUMPTS)	65
Set Minimum Calibration and Validation Values (MIN)	66
Set Maximum Calibration and Validation Values (MAX)	66
Set Oven and Calibrator Delays (DELAY)	66
Set the Calibrator IP Address, Port and Calibrator Number (IPCAL)	67
Set the Calibration and Validation Average (CALAVG)	67
Set Enable CALZ Before validation (VALZO)	67
Set Field Calibration Points (FCAL)	67
Set Field Validation Points (FVAL)	68
Oven Variables (LIST O)	69
List Oven Settings (LIST O)	69
Set the Oven IP Address, Port and Calibrator Number (IPOVEN)	69
Set the Oven Start Command (STARTOVEN)	69
Set the Oven Stop Command (STOPOVEN)	70
Set the Temperature Command (TEMPOVEN)	70
Miscellaneous Variables (LIST M)	71
List Miscellaneous Settings (LIST M)	71
Set Simulation (SIM)	71
Set Echo (ECHO)	71
Set Excitation (XITE)	71
Set Server Select (SVRSEL)	72
FTP Variables (LIST FTP)	73
LIST FTP Settings (LIST FTP)	73
Set User Name (USERFTP)	73
Set FTP Password (PASSFTP)	73
Set FTP Path (PATHFTP)	74
Set FTP IP Address (IPFTP)	74
Set FTP File Name (FILEFTP)	74

SECTION 6: MAINTENANCE	76
Firmware Installation & Configuration File Upload (Windows XP, 7, 10)	76
Calibration Coefficient Installation (Windows XP, 7, 10)	77
DIP Switches & Bootloader Mode	77
DIP Switches	77
When to operate in bootloader mode	77
How to operate in bootloader mode	77
Formatting the Flash Disk	78
APPENDIX	82
Appendix A: Intra-Frame A/D Timing	82
Appendix B - Internal Socket Allocation	83
Appendix C - LabVIEW Client Example	84
Appendix D: 1588 Synchronized Scan Start	85
Appendix E - Engineering Unit Conversion Constants	87
Appendix F - Cable Drawings and Pinouts	88
MPS Premium Power/Serial/Trigger Cable	88
MPS Power/Trigger Cable with Flying Leads	88
MPS Premium Ethernet Cable	89
MPS to RJ45 Ethernet Cable	89
MPS To RJ45 Female	90
MPS4264TCU Cable	90
Appendix G - Software Revision Log	91

[INTENTIONALLY LEFT BLANK]

SECTION 1: SPECIFICATIONS

GENERAL SPECIFICATIONS

Size (WxHxD)		
MPS4264	1.145" x 1.883" x 3.515"	
(including tubes)	[29.08mm x 47.83mm x 89.28m]	
MPS4264TCU	4.295" x 2.350" x 6.882"	
(including tubes)	[109.1mm x 59.69mm x 174.8mm]	
MPS4264TCU with Cooling	4.295" x 2.350" x 7.522"	
(including tubes)	[109.1mm x 59.69mm x 191.1mm]	
Weight		
MPS4264/64CPx	7.65oz [217g]	
MPS4264/64NPx	7.8oz [221g]	
MPS4264TCU	31.52oz [0.89kg]	
Inputs (Px)	MPS4264 0.042" [1.067mm] OD (standard) 0.031" [.787mm] OD (optional) MPSTCU 0.063" [1.600mm]OD (standard)	
Inputs (Cal, Ref, CTL, PRG)	0.063" [1.600mm] OD	
Full Scale Ranges	±4 inH ₂ O, ±8 inH ₂ O, 1psid, 5psid, 15psid [995.4Pa, 1990.7Pa, 6.89kPa, 34.5kPa, 103.4kPa]	
Accuracy		
4 inH ₂ O	±0.20% FS	
8 inH ₂ O	±0.15% FS	
1psid	±0.06% FS	
5psid	±0.06% FS	
15psid	±0.06% FS	
Overpressure (No damage)		
4 inH ₂ O	4psid	
8 inH ₂ O	4psid	
1psid	10psid	
5psid	40psid	
15psid	75psid	
Maximum Reference Pressure	50psig (345kPa)	
Data Output Rate (samples/channel/second)	850Hz (TCP/IP Binary) 2500Hz (Fast Mode)	
A/D Resolution	24bit	

Ethernet Connection 100baseT, MDIX auto-crossing

Power Requirements

MPS4264	9-36Vdc, 3.5W
MPS4264TCU with Heater	24-28Vdc, 45W
External Trigger	5-15Vdc, 6.5mA

Mating Connectors

Ethernet	TE Connectivity PLG 8P8C Mini2
Power	TE Connectivity PLG 8P8C Mini1
MPS4264TCU	

Control Pressure

Requirements	
MPS4264/64CPx	65psi minimum
MPS4264/64NPx	90psi minimum

ENVIRONMENT SPECIFICATIONS

Temperature	
Operation	0 to 70 °C -60°C to 60°C with TCU heater -60°C to 125°C with TCU cooling kit
Storage	0 to 80°C

Humidity	
Operation	5 to 95% RH, Non-Condensing
Storage	5 to 95% RH, Non-Condensing

Shock & Vibration MIL-STD-810G, Category 24

Media Gasses compatible with
Silicon, Silicone,
Aluminum and Buna-N

Maximum Environmental Pressure 100psia (690kPa absolute)

Minimum Environmental Pressure 0.50psia (3.45kPa absolute)

SECTION 1: SPECIFICATIONS

MPS4264

FIGURE 1-1: MPS4264 DIMENSIONAL DRAWING

SECTION 2: INTRODUCTION

OVERVIEW

The MPS4264 is a 64 channel, intelligent, miniature pressure scanner. Each MPS4264 series scanner incorporates 64 individual piezoresistive pressure sensors. In addition to the pressure sensors, the MPS4264 integrates all of the electronics for the analog-to-digital (A/D) conversion process and a processor running a DSP operating system supporting the engineering unit conversion process and all communications overhead. The integral DC/DC converters allow for a single 9-36Vdc power input. Communications and data output is via Ethernet, and RS-232 communications are included for backup and troubleshooting. All 64 pressure sensors share a common reference manifold in the center of the module. The MPS4264 module has two pneumatic valves, each servicing 32 pressure sensors. The valve supports four logic states:

1. Measurement
2. Calibrate
3. Purge
4. Isolate

The valve state is selected by either applying control pressures (applicable to standard units) or via software command (applicable to electric actuation units).

The MPS4264 is the first product designed by Scanivalve to leverage the proprietary “Dynamic CALZ” function. The function effectively nulls the zero offset of the entire analog-to-digital conversion process with every single scan. This feature produces extremely stable sensor outputs and virtually eliminates all zero drift in the system. With negligible zero drift over time or temperature, the traditional pneumatic zero correction that was frequently required is now greatly reduced.

MODULE CONFIGURATIONS

The MPS4264 is the initial pressure scanner in the MPS4000 series of pressure scanners. All variants in the series use the same sensors and valves. However, the electronics are different for each, and the overall form factor is slightly different for each model. Because the interface, form factor and function is different for each variant they will be covered in separate documents. There are several options that, for the most part, apply to all three variants.

PRESSURE RANGE

Each MPS module can be ordered in the following pressure ranges:

4 inH ₂ O	5psi
8 inH ₂ O	15psi
1psi	

VALVE ACTUATION

There are two options for actuating the valves in the MPS4000 series. Each one is optimal for a specific application. They are:

Pneumatic “CPx” - This option places a pneumatic input on each end of the valve. 65psi control pressure applied to each of the ports individually determines the position of the valve. This is optimal for any MPS4000 module that will be replacing legacy ZOC33 modules as the 65psi control pressures are common to both.

Pneumatic “NPx” - This option uses a pneumatic actuator on one end and a spring on the opposing end. The spring is constantly driving the valve into the measurement, or “Px” mode. This requires no control pressure to scan and collect data. The opposing pneumatic control pressure required is 90-120psi.

PRESSURE SENSORS

The 64 pressure sensors in the MPS4000 series are mounted in groups of four to custom designed ceramic bases. Sixteen (16) ceramic bases, or “sticks” are RTV’d to a printed circuit board to create a 32 channel set. This architecture allows for a single four-channel sensor “stick” to be replaced at the factory if needed. It also provides double isolation from any mechanical stresses that may be present in the assembly.

Each sensor PCB includes four (4) digital temperature chips and four (4) RTD's. The RTD's are for use in analog environments (MPS4164) while the digital temperature chips are used in the intelligent environments (MPS4264). The digital temperature chips have been mapped to determine which pressure channels they are most representative of and are used in the pressure engineering unit conversion process.

For every MPS4000 series module, two sensor circuit boards are mounted back to back in the central housing. The core of the housing is an open volume that is used as a common reference pressure for all 64 channels. The design places the sensors deep in the center of the module. Being

located in the center of the module protects the sensors as much as possible from rapid temperature changes and allows for very accurate temperature compensation.

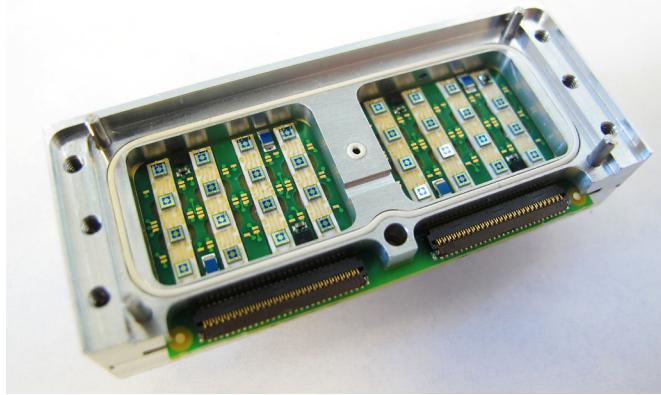


FIGURE 2-1: PRESSURE SENSOR OVERVIEW

PNEUMATIC VALVE

A single MPS4000 series module integrates two pneumatic valves, each supporting 32 channels. The valves are located directly under the pressure inputs on each side of the module. The valve supports four pneumatic states in two physical states. In the “measurement mode” pressures are directed in from the measurement inputs to the individual pressure sensors. The other physical state, “calibrate mode” provides three pneumatic states. They are:

1. Calibrate - A single pressure applied to the “CAL” port is directed to the positive side of all 64 transducers.
2. Purge - Purge pressures applied to the two purge inputs are directed out the 64 measurement ports to clear any debris, particulate or moisture from the input lines.
3. Isolate - The 64 pressure transducers are isolated from the input lines and from the purge pressure.

The valve design is a sliding-type valve. An aluminum shuttle populated with O-rings cycles back and forth between two positions to achieve each of the physical states. An interface between the sliding O-rings and the aluminum stationary portion of the valve is a proprietary compound that is extremely low friction. The O-rings are supported completely to prevent any deformation during a state change. Additionally, the valve shuttle is supported on ball bearings to provide minimal friction and maximum support. This design allows for low actuation forces, minimal “stiction” - meaning the actuation forces do not noticeably increase over time, and minimum wear. Samples of the valve design were tested to over 1,000,000 cycles without any maintenance.

As a subset of the pneumatic option it can either be a

spring-driven “NPx” version or a “CPx” bi-directional pneumatic driven version. The advantage of the “CPx” pneumatic version is that it retains the same 65psi control pressure that legacy ZOC33 analog pressure scanner modules used. Two control pressure inputs are used, Px CTL and CAL CTL. In the case of the spring-driven “NPx” version a spring is constantly driving the valves into measurement, or Px mode. This requires no control pressures to be present for measurement and also prevents the valve from unexpectedly or undesirably moving out of measurement mode. The “NPx” version requires 90-120psi control pressures which is common to the DSA line of products that Scanivalve produces.

All versions of the MPS4264 incorporate an optical sensor used to determine the position of the valve. This can be queried through the MPS’s software interface. The valve position sensor only confirms that the valve is fully in the measurement mode or not. It does not verify with 100% certainty that the valve is in the calibrate mode.

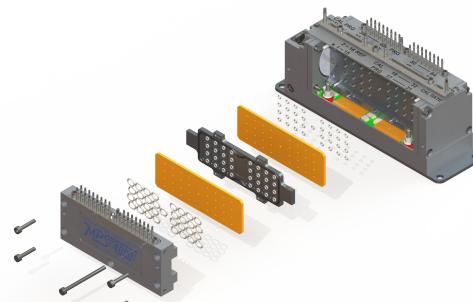


FIGURE 2-2: PNEUMATIC VALVE EXPLODED VIEW

ELECTRONICS

The MPS4264 has seven (7) unique printed circuit boards including two (2) flexible circuits. Under the cover on the end of the module is the IO Board (Input-Output) which supports the power input and communications. This board has all of the interface electrical connectors including the Ethernet, power and RS-232 connectors. It is also where the DC-DC power supply is located. The IO Board is connected to the processor board via a flexible circuit. The Processor Board is where the modules memory and processor reside. This board is located along the very bottom of the module inside of the lower aluminum cover. Above the Processor Board is the A/D board which is connected by a rigid parallel board-to-board connector. The A/D board is where all 16 A/D converters are located (one A/D per 4 transducers) as well as all of the multiplexing and sensor excitation. The A/D board is connected to each of the two Sensor board with two flexible circuits each. Each connection handles a quadrant of 16 channels. The Sensor board, as discussed earlier has 32 individual pressure sensors mounted on ceramic bases. The last circuit board is the valve position sensor board. The single valve position sensor resides in a machined cavity of the valve/side plate and optically checks the position of the valve shuttle. It connects directly to the processor board using a miniature connector.

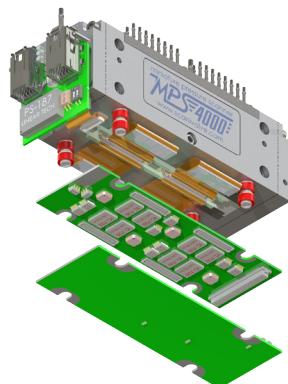


FIGURE 2-3: ELECTRONICS EXPLODED VIEW

PNEUMATIC INPUTS

For added convenience and flexibility, all pneumatic inputs on the MPS4000 series use removable headers. There are three separate manifolds, two headers for the 64 input ports and a third manifold for the CAL, REF and Control Input(s). The purge supply pressure is also applied through each of the two measurement manifolds. The manifolds are held on with three captive screws each. Under the CAL/REF/CTL manifold is an ID plate that will be free when you remove the manifold. Take care to not lose this plate or the O-rings that are retained by it.

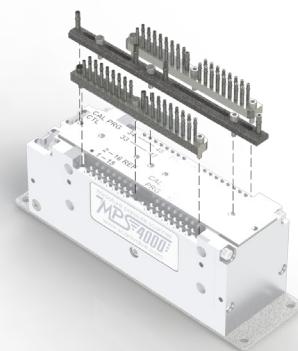


FIGURE 2-4: PNEUMATIC INPUT HEADERS

POWER & ETHERNET CONNECTIONS

The MPS4264 has its power and Ethernet connections located on one end. The power connection also serves as an external trigger input and a serial connection. LEDs are positioned in a recessed cavity between the two connectors to indicate power, link and activity. Both the power and the Ethernet connector are TE Connectivity "Mini IO" series but are polarized differently to prevent mistakes. The connectors are latching to prevent any vibration from loosening the connections.

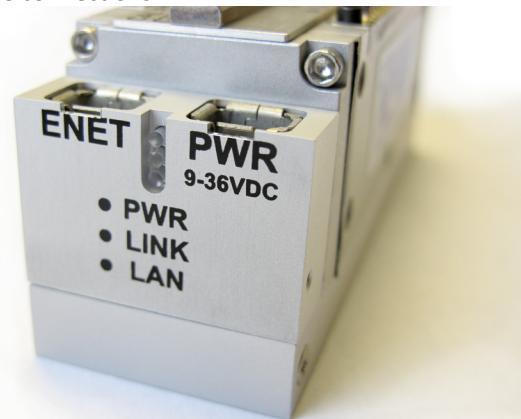


FIGURE 2-5: POWER & ETHERNET CONNECTIONS

[INTENTIONALLY LEFT BLANK]

SECTION 3: OPERATION

UNPACK & INVENTORY

When you first unpack the MPS4264 module, begin by inspecting and inventorying the contents of the package. If any visible damage is immediately noticed or if any contents are missing, contact Scanivalve before proceeding. As a minimum, MPS4264 modules are shipped with the following contents:

1. MPS4264 module
2. Protective case
3. 155625 Power cable with flying leads
4. 156110 Ethernet RJ45 adaptor cable
5. MPS4264 resource CD

MOUNTING

The MPS4264 includes a stainless steel mounting plate. Using three #2-56 screws it can either be mounted to the bottom of the module or to the side. The mounting plate is 0.060" (1.5mm) thick with four 0.125" (3.1mm) holes in it. See Figure 3-1 for hole spacing, and Figures 3-2 and 3-3 for instructions for attaching the mounting plate to the MPS module. The MPS can be mounted in any orientation.

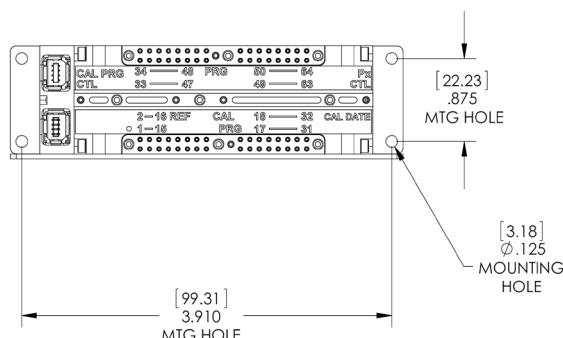


FIGURE 3-1: MOUNTING PLATE DIMENSIONS

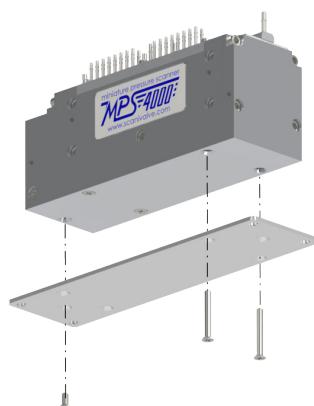


FIGURE 3-2: MOUNTING PLATE, BOTTOM

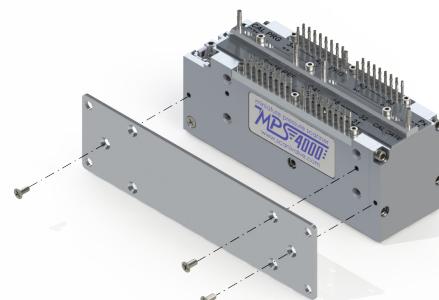


FIGURE 3-3: MOUNTING PLATE, SIDE

ENVIRONMENTAL CONSIDERATIONS

The MPS4264 series is designed to withstand normal industrial, flight test, educational, wind tunnel or similar applications. The module is not water proof and it should be protected from any splash, spray or mist. If any moisture gets spilled or splashed on the MPS module, wipe it dry immediately to prevent damage to the module.

The MPS module should not be mounted in a location where it may be subjected to extreme temperature shifts or ambient temperatures outside limits defined in "Section 1: Specifications" on page 7. Keep in mind that the internal temperature of the module will run approximately 15°C - 20°C warmer than ambient temperature but this is accounted for in the calibration tables. The temperature specification is with regards to the ambient temperature.

If the MPS module needs to be used in an environment outside of the guidelines above it should be installed in a Thermal Control Unit (TCU). This is a rugged enclosure incorporating a heater or cooling fittings in order to protect the MPS module.



CAUTION! Mounting the MPS module inadequately or in an environment that does not conform to the recommendations can result in permanent damage to the module.

WARM-UP

Because of the design of the MPS4264, the accuracy is only minimally impacted by the device warming up after initially being powered. Figure 3-4 below shows a sample 5psi module being powered on. Pressures from 3 channels are shown along with the module temperature. As indicated in the graph, the first 15 minutes of warm up produces the only noticeable change in the pressure readings. The module can take up to 3 hours to fully warm up (in a 25°C ambient environment) but throughout that warm-up the pressure readings are very good. It is recommended that the module be allowed to warm up for a 15-30 minutes before collecting data, but this is not mandatory. If time allows, the warm-up period should be extended to 60 minutes for most applications.

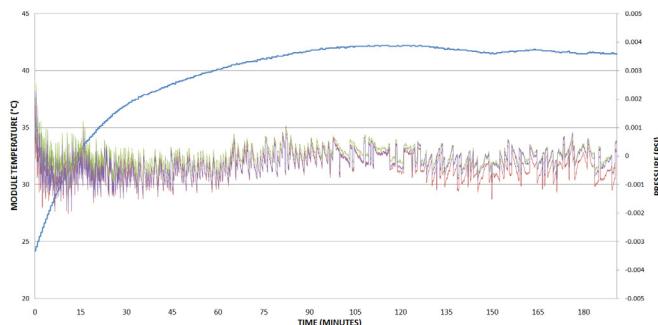


FIGURE 3-4: TEMPERATURE & PRESSURE THROUGH MODULE WARM-UP

COMMUNICATIONS

The MPS module is designed primarily for Ethernet communications. This provides a means to configure the MPS module as well as scan and collect data from the module.

A serial RS-232 connection is also supported. The serial connection is designed to provide emergency communications with the module in cases where an Ethernet connection cannot be established. If the internal memory in the MPS is corrupted in any way, the serial connection allows the user to format the memory and restore the operating files.

SERIAL COMMUNICATIONS

Every MPS4264 module has an RS-232 serial output. It is available through power connector on pins 1, 2 and 3. The wiring diagrams can be found in “Appendix F - Cable Drawings and Pinouts” on page 88. Settings for establishing a serial connection to the MPS4264 module are as follows:

Bits per second: 9600 BAUD

Data bits: 8

Parity: none

Stop bits: 1

Flow control: none

Once a connection is established any of the commands listed in “Section 5: Software” on page 33 can be sent to the module. The only limitation is that a SCAN cannot be initiated via the serial connection and scan data will never be returned via the serial connection. Most importantly, through the serial connection, the module’s “Group IP” settings, including the IP address, can be viewed and set. The module’s internal memory can also be formatted and restored through a serial connection (see “Formatting the Flash Disk” on page 78 for this procedure.)

Two power/serial cables are available for the MPS. They are:

P/N 156085: “Premium” Mini I/O male to Mini I/O male power/serial/trigger cable

P/N 155625: Mini I/O Trigger/Power cable with flying leads

The cable listed as “Premium” versions use a highly-flexible, silicone coated cable that is the smallest possible diameter (0.156”, 4mm).

ETHERNET COMMUNICATIONS

The MPS4264 has one Ethernet connection, 10/100Base-T with MDIX auto-crossing. This is the primary means of communications with the MPS. All MPS4264 modules use a TE Connectivity “Mini I/O” series connector for the Ethernet connection. This is a latching connector to prevent the connector from becoming unseated when subjected to vibrations. Shielded Category 5e cable or better is recommended for all Ethernet connections. Several variations of Ethernet cables using the TE Connectivity “Mini I/O” connector are available from Scanivalve. Depending on the application, different cables may be required.

ETHERNET CABLES:

P/N 155635 Mini I/O male to RJ45 male patch cable

P/N 156062: “Premium” Mini I/O male to RJ45 male patch cable.

P/N 156110: Mini I/O male to RJ45 female receptacle extender cable.

P/N 156091: “Premium” Mini I/O male to Mini I/O male patch cable (MPS - ES4000 Switch)

Cables listed as “Premium” versions use a highly-flexible, silicone coated cable that is the smallest possible diameter (0.156”, 4mm). Maximum cable length for any of the Ethernet cables is 200 feet. For additional information on establishing an Ethernet connection with the MPS, see “Setting Up an Ethernet Connection” on page 1616.

CLIENT/HOST OPTIONS

Once the module has been connected and the IP address has been configured, communications can be established with the MPS module. Communications can be made through several software packages including:

- PC - TCP/IP
- PC - UDP
- PC - ScanTel (Scanivalve PN: 155406-01)
- PC - LabVIEW
- PC - Windows HyperTerminal

PC - TCP/IP

The user may write their own TCP/IP interface using the software specification portion of this manual. This interface should allow the user to:

- Issue commands to any or all MPS modules on the network.
- Display returned information or scan data from the MPS module(s).
- Write returned information or scan data to the client/host in TCP/IP format.

PC - UDP

The MPS provides a dual-output function in which data can be output simultaneously out via a FTP connection and a UDP broadcast data stream. This allows for great flexibility in the configuration of the data collection. It allows users to configure a system such that data can be monitored “real-time” via the UDP stream and collected reliably via the FTP data transfer.

PC - SCANTEL

ScanTel is a free communications utility designed by Scanivalve to communicate with Scanivalve products including MPS4264 modules. It is a text based, command line program that allows users to connect to a single MPS module and modify the configuration variables, upload or download coefficients and collect data in both TCP/IP and UDP format.

PC - LABVIEW

The MPS4264 has been developed specifically with a simplified LabVIEW interface in mind. This allows the user to easily create a custom graphical interface with the MPS to configure the unit, collect and display data. For detailed information on setting up a LabVIEW interface with the MPS, see “Appendix C - LabVIEW Client Example” on page 84. Free sample VI’s are available to assist with creating a custom LabVIEW interface.

PC - HYPERTERMINAL

HyperTerminal is a Windows program included as part of Windows 2000, XP and Vista Operating Systems. This program permits a user to connect to a single MPS module, modify the configuration variables, upload or download

coefficients and collect data. HyperTerminal provides a means for both Serial RS-232 and Ethernet connections. It is a text based command line program. HyperTerminal emulators like “PuTTY” can be downloaded for free from various sources online if HyperTerminal is not available.

FTP

The MPS4264 can transfer scan data via FTP to a FTP server. The user may set up any suitable server or NAS capable of running as a FTP server. This provides a secure, reliable data connection to generate data files directly on an FTP server.

SETTING UP AN ETHERNET CONNECTION

Before an Ethernet connection can be established the IP address need to be configured. This can be done manually with a “static” IP address or automatically using a “dynamic” IP address automatically assigned to the MPS by a DHCP server. If you will be using a “static” IP address, then the IP address of the MPS must be manually set by the user. The IP address assigned to the MPS must be compatible with the network/host computer. If you will be using a “dynamic” IP address, the MPS must be both configured to obtain an IP address automatically and it must be on a network with a DHCP server present.

USING A STATIC IP ADDRESS

If you choose to use a static IP address, you must manually set the IP address of the MPS. The IP address you choose to use must be compatible with the network the MPS is being connected to. If you are simply connecting the MPS directly to a host computer, then the IP address of the MPS and the host computer must be compatible and should both be configured manually. The range of compatible IP addresses is defined by the subnet mask. The standard default subnet mask is 255.255.0.0. This default subnet mask requires the IP address of the module and host computer must share the first two octets, or sets of numbers. The third and fourth octets of the IP address is variable, although it is recommended that the third octet also be shared between the host computer and the module. The subnet mask digits of “255” define that the two IP addresses must have matching digits in those positions and the subnet mask digits of “0” allow the two IP addresses to have unique values for those octets. No two devices on a single network can share the same IP address. Below are some examples of compatible IP addresses.

Example of matching the first three octets:

Subnet mask: 255.255.255.0
Host computer: 191.30.90.100
MPS module: 191.30.90.125

Example of matching the first two octets:

Subnet mask: 255.255.0.0
Host computer: 191.30.1.100
MPS module: 191.30.90.125

Example of NON-COMPATIBLE IP addresses:

Subnet mask: 255.255.255.0
Host computer: 191.30.5.1
MPS module: 191.30.90.5

When connecting the MPS and configuring the network, you can choose to either change the IP address of the MPS module to match the network, or change the IP address of your host computer to be compatible with the MPS module. If you choose to change the IP address of the host

computer to be compatible with the MPS, then:

The IP address of a Windows XP host computer can be changed under:

Control Panel -> Network Connections -> Local Area Network -> Properties -> Internet Protocol (TCP/IP) -> Properties.

The IP address of a Windows 7 host computer can be changed under:

Control Panel -> Network and Sharing Center -> Local Area Connection -> Properties -> Internet Protocol Version 4 (TCP/IPv4) -> Properties.

MPS4264 modules are shipped with a preset IP address that uses the following format:

MPS4264: 191.30.90.XXX

(Where “XXX” is the last three digits of the unit’s serial number.)

The MPS4264’s IP address is set with the “IPADD” variable in the IP Group of variables. This can be changed through an existing Ethernet connection or through a serial connection. If an existing Ethernet connection can be established (by changing the host computer’s IP address or some other means) it is the preferred method. If it is not possible, then a serial connection can be established with the MPS. **If the IPADD variable is changed either through an Ethernet connection or a serial connection, the new address will not be effective until a “SAVE IP” command has been issued and power has been cycled.** For more information on the IPADD variable, “Set IP Address of the MPS (IPADD)” on page 59. More information on establishing a serial connection can be found in “Serial Communications” on page 14. Once the IP address of the MPS and host computer have been set to be compatible, the connection should be tested. The best means of doing this is to “ping” the MPS using the host computer’s command prompt.

USING A DYNAMIC IP ADDRESS

The MPS can be setup to have a IP address dynamically assigned from a DNS server using DHCP. When the IP address is set to 0.0.0.0 then DHCP is used. When DHCP is active, the MPS’s DNS name is used to connect to the unit. The DNS name is defined by the variable NAME in the LIST IP group.



Remember that any changes to the variables in the LIST IP group require a SAVE IP command be issued and the MPS must be restarted for them to take effect.

The DHCP server will supply an IP address to the MPS for a period of time, called a lease. The MPS continues to broadcast discovery packets every 5 seconds until the DHCP server responds with an “offer.” Network communication is not possible to the MPS unless it is in the “IP is Now Leased” state. See Figure 3-6 below. The lease is usually granted for 24 hours. This sequence also applies to renewals that typically occur every 12 hours.

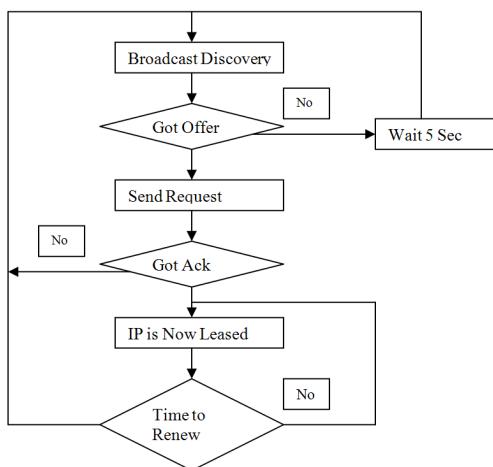


FIGURE 3-7: DHCP LEASING PROCESS

The MPS only issues a discovery packet to the DHCP server. It does not communicate with any DNS server that may be attached to the network. Therefore when it is desired to use a name in place of an IP address, the network must be set up so that the DHCP server retrieves the name from the DNS server. This is referred to as DNS dynamic record update. The MPS does furnish a requested name to the DHCP server so that it may pass the name on to the DNS server.

For Windows Server based networks, Figure 3-8 shows the settings for automatic dynamic update of the DNS server. The exact method for making these setting will vary based on the operating system and operating system version.

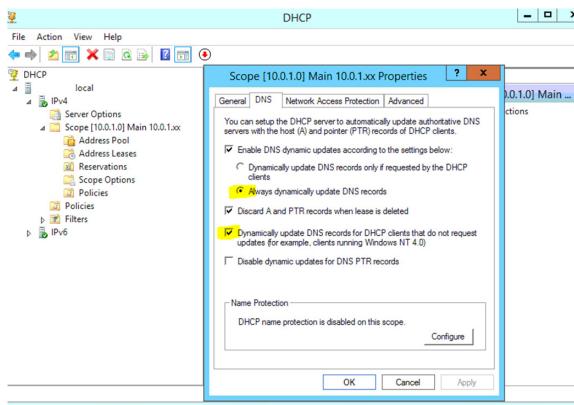


FIGURE 3-8: WINDOWS SETTINGS FOR AUTOMATIC DNS SERVER UPDATES

Once the IP address of the MPS and host computer have been set to be compatible, the connection should be tested. The best means of doing this is to “ping” the MPS using the host computer’s command prompt.

WEB SERVER: CONNECTING

The MPS4264 is a single user web server which provides a simple method for communication, configuration and data collection from virtually any host computer. It requires no additional software be installed on the host computer and provides a simple graphical interface for the user.



FIGURE 3-9: WEB SERVER WELCOME SCREEN

The web server has been tested with several web browsers to verify functionality and compatibility. Those that were tested are:

Internet Explorer 10 (and later)
Microsoft Edge

To connect to the MPS4264’s web server, simply open one of the supported browsers. In the address bar, enter:

xxx.xxx.xxx.xxx/index.htm
(where xxx.xxx.xxx.xxx is the IP address of the MPS).

Once connected, the web server’s welcome screen will be displayed.

At the top right corner of the screen, basic module information will be displayed. This will include IP address, Model Number, Serial Number and Nominal Pressure Range.

WEB SERVER: CONFIGURING SCAN SETTINGS & SCANNING DATA TO A FILE

Once you have connected to the MPS4000 and have the home page displayed, you may click on various links to perform different tasks. On the blue navigation bar select the “Scan” link. You will get the following screen:

SECTION 3: OPERATION

MPS4264

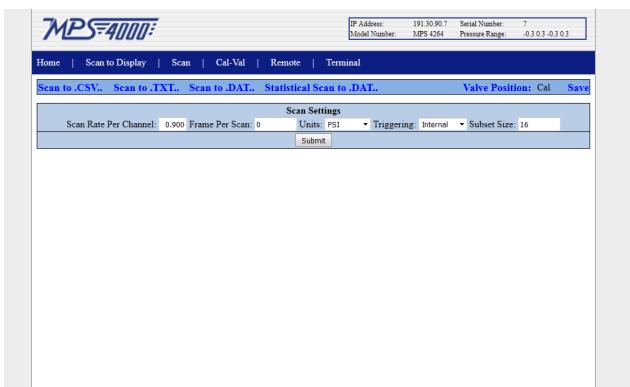


FIGURE 3-10: WEB SERVER SCAN CONFIGURATION SCREEN

There are five configuration variables that you can set from this screen. They are:

- SCAN RATE PER CHANNEL: This sets the data output rate in samples/channel/second (Hz).
- FRAME PER SCAN: This sets the number of samples that will be output for each time a “SCAN” command is received.
- UNITS: This is a selection box, allowing you to select the pressure units the data will be output in.
- TRIGGERING: This is a selection box, allowing you to select the method that a each frame of data is initiated by. There are four options:
 - INTERNAL - A scan is initiated, and subsequent frames of data are scheduled using the MPS’s internal clock.
 - EXTERNAL - A scan is initiated, and subsequent frames of data are schedule using an external trigger.
 - AUTOSCAN - The MPS will begin scanning and outputting data immediately on power up. See “Auto Scan Option” on page 24 for more information.

Additionally, along the top of the screen there are four “Scan to...” hyperlinks allowing you to initiate data collections using different methods. Those four options are:

- SCAN TO .CSV: This generates a .CSV file that will be saved to your local computer. This file will contain ASCII data in a column format, with a column of data for each pressure and temperature channel.
- SCAN TO .TXT: This generates a .txt file that will be saved to your local computer. This file will contain ASCII data in a single “scrolling” format, with all data from each frame of data printed out below the previous.
- SCAN TO .DAT: This generates a .dat file that will be saved to your computer. This file will be binary data that will need to be converted to ASCII as a post processing task.
- STATISTICAL SCAN TO .DAT: This generates a .dat file

including pressure data and statistical analysis data that will be saved to your computer. This file will be binary data that will need to be converted to ASCII as a post processing task. See “Statistical Scan Option” on page 24 for more information.

On the same header bar as the “Scan To...” hyperlinks are two additional options. The “VALVE POSITION:” is an indicator identifying the state of the pneumatic valve in the MPS. It will either be “Cal” indicating the module is in calibrate mode or “Px” indicating the module is in measurement mode. **This is not automatically updated when the valve state change, it must be manually polled. Do this by single clicking on the blue “Valve Position” hyperlink.** The final option in the header bar is the “SAVE” option. Clicking this issues a SAVE command to the MPS taking the current configuration from volatile RAM memory and saving it in the non-volatile FLASH memory storage.

At the bottom of the screen is the “SUBMIT” button. For any changes that are made to be sent to the MPS and take effect, the “SUBMIT” button must be clicked. Without clicking “SUBMIT” after making a change, the command will not be sent to the MPS and the change will not be effective.

WEB SERVER: SCANNING DATA TO THE SCREEN

After you have set up the scanning configuration you are now ready to scan data to the display. On the blue navigation bar, click on the “Scan To Display” link. You will get the following screen:

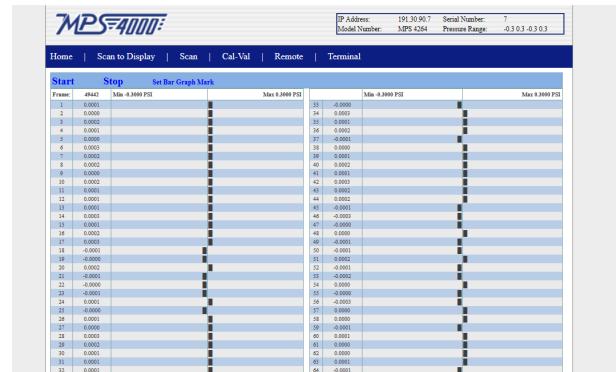


FIGURE 3-11: WEB SERVER SCAN TO DISPLAY

Click on the “START” link. The bar graphs will display the pressure as frames of scan data are sampled. Click on “STOP” when finished. There are three important notes about the “Scan to Display” function and how the bar-graph works.

- 1) The bar graph will never update faster than 1 Hz. If the MPS is set to scan at a higher rate than this, the data collected between the display refreshes are ignored and discarded.
- 2) When the “Scan to Display” function is used, the “FPS” variable is ignored. The MPS will continue to scan and the bar graph will display data until the STOP button is clicked no matter what the FPS variable is set to.
- 3) The bar-graph automatically scales based on the highest pressure value that is being displayed.

The size and scale of the “marks” in the bar graph can be set by clicking “Set Bar Graph Mark.” This only effects how the bar is displayed in the bar graph and has no effect on the data or scaling of the graph.

WEB SERVER: CALIBRATIONS & VALIDATIONS

THIS FEATURE WILL BE RELEASED IN A FUTURE SOFTWARE UPDATE. Calibrations and validations of the MPS can be controlled through the web server as well. The web server’s calibration-validation interface makes setting up a calibration and/or a validation much simpler than manually setting up each of the individual parameters for the operation. It is accessed by clicking on the “Cal-Val” link in the blue navigation bar.

Calibration Temperature	Validation Temperature	Calibration Pressure	Validation Pressure
Temperature points: 15 Minimum temperature point: 0.0000 Maximum temperature point: 79.0000 Temperature block size: 191.30.40	Temperature points: 3 Minimum temperature point: 10.0000 Maximum temperature point: 60.0000	Pressure points: 25 Minimum pressure point: -0.3179 Maximum pressure point: 0.3179 Pressure block size: 90	Pressure points: 15 Minimum pressure point: -0.2890 Maximum pressure point: 0.2890
Oven IP: _____ Port: _____ Dwell (Minutes): _____ Calibrator IP: _____ Port: _____ Number: _____ Dwell (Seconds): _____			

FIGURE 3-12: WEB SERVER CAL-VAL SCREEN

The Cal-Val screen is broken up into three primary sections: Temperatures, Pressures and Communications. The “Temperatures” section includes several variables, each duplicated to allow independent configurations for Calibrations and Validations:

TEMPERATURE POINTS: This is used to set the number of temperatures used during the calibration or validation. If the “temperature points” is set to 0, the MPS will do a calibration/validation only at the module’s current temperature and will not control the oven.

MINIMUM TEMPERATURE POINT: This is used to set the lowest temperature point in the calibration/validation

MAXIMUM TEMPERATURE POINT: This is used to set the maximum temperature point in the calibration/validation.

The MPS will use the minimum temperature and maximum temperature and divide the spread evenly using the number of temperature points to determine the temperatures that the calibration/validation will be performed at.

The “Pressures” section includes several variables, each duplicated to allow independent configurations for Calibrations and Validations:

PRESSURE POINTS: This is used to set the number of pressures that will be applied during the calibration or validation.

MINIMUM PRESSURE POINT: This is used to set the lowest pressure applied during the calibration/validation.

MAXIMUM PRESSURE POINT: This is used to set the maximum pressure applied during the calibration/validation.

The MPS will use the minimum pressure and maximum pressure and divide the spread evenly using the number of pressure points to determine the pressures that are applied to the module during the calibration/validation.

The “Communications” section includes parameters for establishing a connection between the MPS, the pressure calibrator and the oven (if applicable):

IP (OVEN): Sets the IP address of the oven controller.

PORT (OVEN): Sets the port number of the oven controller.

DWELL (MINUTES) (OVEN): Sets the dwell time, in minutes to allow the oven to stabilize at a temperature after setting the temperature.

IP (CALIBRATOR): Sets the IP address of the pressure calibrator.

PORT (CALIBRATOR): Sets the port number of the pressure calibrator.

NUMBER (CALIBRATOR): Sets the channel number of the calibrator being used.

DWELL (SECONDS) (CALIBRATOR): Sets the dwell time, in seconds to allow the calibrator to stabilize at a pressure. If this is set to 0, the MPS will query the calibrator to determine stability and advance as soon as the calibrator is stable.

Across the top of the screen is a bar with several links/options. These features allow you to start and stop a calibration, as well as interface with the module to retrieve the latest validation report, download a copy of the calibration

coefficients and perform some checks. These options are:

START: Click this to initiate a calibration/validation sequence. This sequence will include a calibration and/or a validation as defined by the options that are enabled with the “Calibration” and “Validation” check boxes.

STOP: Stop the current operation, whether a calibration or a validation.

GET REPORT...: Download the latest validation report.

This will be saved as a .CSV file called “Report.CSV”

GET TABLE...: Download the current calibration coefficients. This will be saved as a .txt file called “Insert.TXT”

CAL/VAL POINTS CHECK: Performs a check and outputs the temperatures and pressures to be used during the calibration/validation based on the pressure and temperature min/max variable settings.

DELETE ALL: Deletes the current calibration coefficients from the module.

ZERO CAL: Performs a zero correction calibration, called “CALZ.” Any pressures applied to the module during the CALZ will be subtracted out as offset.

SAVE: Issues a “save” to the module saving the current settings and calibration coefficients.

Once all of the Calibration/Validations configurations have been set, click the ‘Submit’ button at the bottom of the Cal-Val screen. At this point, clicking ‘Start’ in the header bar will initiate the calibration and/or validation sequence. If a calibration is being performed, then new calibration coefficients will be generated. These will be uploaded to the module but may or may not over-write the existing calibration coefficients depending on the exact setup and the previous calibration coefficients. Before beginning a calibration, the old coefficients should be deleted. This is done by clicking the ‘Delete All’ link in the header bar. If it is desired, the existing coefficients can be backed up before being deleted using the ‘Get Table...’ link. It is recommended that after a calibration has been performed, the new calibration coefficients should be retrieved from the MPS and backed up in a safe location.

After a validation has been performed, the validation report will be generated and saved to the flash memory of the MPS4264. This report can be retrieved using the ‘Get Report...’ link. The report will be downloaded as a .CSV file for easy viewing in Microsoft Excel or similar.

ZERO CAL

The MPS4264 was designed for maximum sensor stability and extensively tested for operations with no zero correction calibrations. For many applications the ‘Zero Cal’ (aka CALZ) feature will not be needed. Every sensor behaves in a unique way so it is impossible to outline specific guidelines on when the CALZ will need to be used, but for most cases where the module is operating between 15°C and 40°C no CALZ should be needed. Scanivalve recommends monitoring and checking the module’s readings with zero pressure applied to determine if sensors have drifted and if a CALZ needs to be performed.

WEB SERVER: “REMOTE” SCREEN

The MPS4264 incorporates a feature to allow multiple MPS4264 units to communicate with one another, allowing certain commands to be sent to a single “master” unit and then rebroadcast out to the selected “slave” units. This makes the user interface much simpler and allows for more accurate synchronization between units. This function uses the network multicast feature to establish an MPS4264 group or “cluster.” While the multiple-unit commands and interface is limited through the web server, the “Remote” tab on the web server allows you to search for any MPS4264 scanners that are in the same “cluster” as the MPS you are currently communicating with is.

To search for MPS scanners that are in the same “cluster” as your master device, simply click “Find Remote Devices” in the upper navigation bar. A pop-up screen will appear advising you that the master is searching for slave units. Once the pop-up closes, you will be returned to the web server. If any modules are discovered in the cluster, they will be listed including their serial number, IP address and pressure range as shown in Figure 3-13 below.

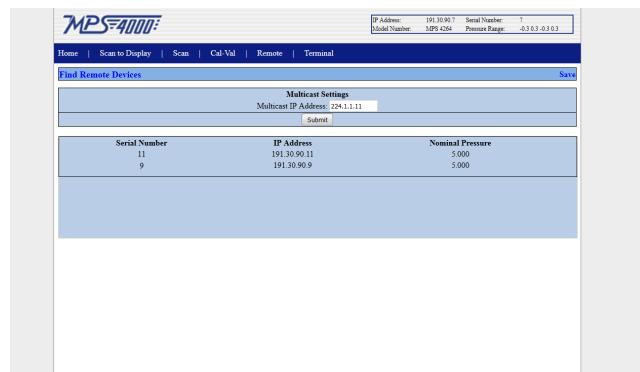


FIGURE 3-13: WEB SERVER REMOTE SCREEN

In addition to searching for other MPS4264’s in the cluster, the “Remote” page of the web server allows you to change the unit’s current multicast address. This can be set in the

“Multicast IP Address” field. In order for the change to take effect, you must:

1. Enter the desired multicast address in the field
2. Click the “submit” button
3. Click the “Save” link at the top of the screen
4. Wait for the save to complete (typically 30 seconds)
5. Cycle power on the MPS.

Once this has been completed the MPS4264 will be configured with the new multicast address.

For more information on the multicast function and working with multiple MPS units at the same time, see “Multiple Unit Operations” on page 22.

WEB SERVER: “TERMINAL” SCREEN

The MPS4264’s web server integrates a “terminal” emulator screen allowing the user to send simple text commands to the scanner. While most of the configuration variables in the module are accessible through the several GUI web server screens, not every variable in the MPS can be queried or changed through the graphical portions of the web sever. For the remainder of the variables that are not accessible, the “terminal” screen provides a simple means to check and change every variable in the MPS. LIST commands can be sent to check the current setting of relevant variables and any command can be sent to the MPS through the web server’s terminal.

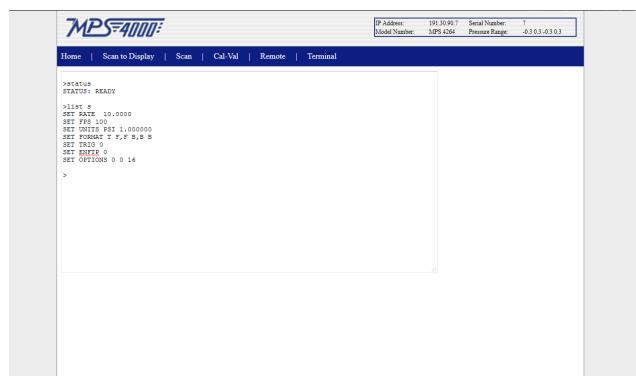


FIGURE 3-14: WEB SERVER TERMINAL SCREEN

The only limitation to the Terminal is that a SCAN cannot be initiated through the Terminal. Because the SCAN command cannot be sent, data cannot be displayed to the terminal or collected using the terminal.

FTP SCAN DATA TRANSFER

The MPS4264 can transfer scan data via FTP to a FTP server. The user may set up any suitable server or NAS capable of running as a FTP server.

See the section on FTP settings to configure the MPS4264 for FTP transfer. The conditions of SSD and SST must be satisfied before scanning will start.

When ENFTP is set to one, and a scan is started, the software will attempt to open a file on the FTP server as configured by all of the FTP settings. The time and date that is indicated by the GetTime command is used in the construction of the file name.

The file that is saved on the server is assembled as follows:

<path on server><base data filename><date>_<time>.<type>

Where <type> is .DAT - binary, .TXT - ASCII Text, .CSV - ASCII CSV file.

The file is closed when the scan stops. If the file cannot be opened on the server when the scan starts, the scan is stopped before any data is output.

When collecting binary data to a FTP server with a scan rate above 200Hz, the MPS will buffer 4 frames of data before transmitting the data. In order to receive all frames in a scan, FPS will need to be set to a value that is divisible by 4. You will not receive the last frames of a scan (maximum total loss of 3 frames of data) if the total amount of frames collected is not divisible by 4.

Example: if FPS is set to 102, and the scan rate is 500Hz, you will receive frames 1 - 100, but you will not receive frames 101 and 102.

THE BINARY TCP SERVER

The MPS4264 provides a binary TCP server designed to support LabVIEW or any other computer based client. The MPS expects a connection to arrive on port 503. **ONLY ONE CONNECTION IS SUPPORTED FROM THE SERVER** due to socket limitations in the MPS. Additional attempts to connect to port 503 when an existing connection is present will be refused.

To enable the binary server the SVRSEL must be set to 2. When the SVRSEL is changed, the group must be saved with a SAVE M command and the MPS rebooted.

The format of the binary packet sent from the binary

server is based on the setting of the FORMAT configuration variable parameter B. See “Set Scanning Data Format (FORMAT)” on page 55 description for more information.

COMMANDS

Because all of the configuration settings may be handles via the Telnet server (port 23), the binary server only supports starting and stopping a scan. Sending a 1 will start the scanning and sending a 0 will stop the scanning. The scanning will also stop if the MPS scan buffer overflows or when the binary client disconnects from the MPS.

Additionally, when a client is connected to the binary server port, a scan command issued via Telnet port causes scanning to start and data to be sent out the binary port.

SCAN DATA

Data sent over the binary server can have three formats: LabVIEW, binary packet or statistical binary packet. Binary packet is described in “Table 5-5: Binary Data Format” on page 37 and is configure via SET FORMAT B B. Statistical binary packet is described in “Table 5-6: Statistical Binary Data Format” on page 38 and is configured via SET FORMAT B S. The following section describes the SET FORMAT B L configuration.

LABVIEW SCAN DATA

When FORMAT B is set to L, scan data, arranged for LabVIEW, is sent from the MPS to the client in binary form that includes the frame number, the MPS average temperature and the 64 channels of pressure data. All values are sent in floating point format making it easier to handle by the client. Each frame contains 264 bytes or 66 array elements. The table below shows each frame’s format:

Parameter	Number of bytes	Array element
Frame Number	4 bytes formatted as float	0
MPS Average Temp	4 bytes formatted as float	1
Pressure Channel 1	4 bytes formatted as float	2
Pressure Channel 2	4 bytes formatted as float	3
.....
Pressure Channel 64	4 bytes formatted as float	65

LABVIEW DATA BUFFERING

The binary server can buffer up to 170 frames of data before the buffer overflows and stops scanning. This allows the host software to gather data from the MPS at a slower rate than the MPS is scanning. However, multiple frames must be read, at once, with each read.

To calculate the maximum time that the client can wait before reading data from the binary server use the following:

$$\text{Max Delay} = (1/\text{MPS Rate}) * (\text{Max Frames in Buffer}) / 2$$

$$\text{Max Delay} = (1/\text{MPS Rate}) * 170 / 2$$

$$\text{Max Delay} = 85 / \text{MPS Rate}$$

Note: The buffer size in the MPS 4000 is 170 frames. However, to handle delays in the network connection only half should be used as a target delay. Thus 85 frames are the maximum frames to hold in the buffer.

EXAMPLE LABVIEW

“Appendix C - LabVIEW Client Example” on page 84 shows an example LabVIEW client interface to the MPS4264 binary server.

MULTIPLE UNIT OPERATIONS

When working with multiple scanners, the operation can be greatly simplified by using the integrated “multiple unit” architecture. When the desired units are configured into “clusters” of modules, a single connection can be established to any one of the modules in the cluster. Once this connection is established, any of the “Multiple Unit” commands can be sent to the master unit which will in turn pass of the command to the other modules in the cluster. This allows many unit to be started scanning, stopped or calibrated all at the same time from a single command sent to a single MPS4264 unit.

The MPS4264 uses the network multicast feature to establish a MPS4264 device cluster. In all cases the user connects to one MPS4264. This device becomes the master device. This master then passes the relevant commands on to the other devices that are members of the device cluster. Each MPS4264 has the ability to be a master or slave.

The membership of a MPS4264 into a cluster is established by setting the variable MCAST to the same multicast address. All MPS4264 devices with the same MCAST setting become members of the same MPS4264 cluster. Note that changing the MCAST variable requires a SAVE to be issued to the unit and the MPS to be rebooted before the change becomes effective. Multiple clusters may exist on the same network. Multicast addresses are limited to the range from 224.0.0.0 through 239.255.255.255.

Once the MPS4264’s have been configured to be included in a cluster and a connection has been established with the “master” unit - issue the command:

MFIND

A list of MPS modules included in the cluster will be

returned. Review the list to ensure it is inclusive of all units that are expected to be included in the cluster. Once this has been verified, a scan on all units can be initiated with the command:

MSCAN

If at any point the scan needs to be terminated, issue the command:

MSTOP

There are three “multiple unit” commands that allow control of the cluster of scanners. They are:

MULTIPLE UNIT SCAN COMMAND

Multiple MPS units are started scanning by connecting to one MPS and issuing a MSCAN command. This device becomes the master and reissues a multicast SCAN command. See the table below for the scan start packet format. All devices that are set up to listen to this multicast address will respond to the SCAN command. The slave devices must still fulfill the SSD and SST requirements for starting a scan.

TABLE 3-1: MSCAN PACKET FORMAT

FUNCTION	BYTES	DATA TYPE	DESCRIPTION
Packet Type	2	Integer	0202

Further information on the MSCAN command can be found in “Multiple Unit Scan” on page 50.

MULTIPLE UNIT STOP COMMAND

Multiple MPS units are stopped scanning by connecting to one MPS and issuing a MSTOP command. This device becomes the master and reissues a multicast STOP command. See the table below for the stop packet format. All devices that are set up to listen to this multicast address will respond to the STOP command.

TABLE 3-2: MSTOP PACKET FORMAT

FUNCTION	BYTES	DATA TYPE	DESCRIPTION
Packet Type	2	Integer	0303

Further information on the MSTOP command can be found in “Multiple Unit Stop” on page 42.

MULTIPLE UNIT FIND COMMAND

The user may connect to one MPS4000 and issuing a MFIND command. This command will return a list of all devices that are part of the multicast cluster. The MFIND command used the Discovery Packet described below.

TABLE 3-3: MFIND PACKET FORMAT

FUNCTION	BYTES	DATA TYPE	DESCRIPTION
Packet Type	2	Integer	0101

Further information on the MFIND command can be found in “Multiple Unit Find” on page 42.

SCANNING WITH AN EXTERNAL TRIGGER

The MPS4264 can be configured to scan with an external hardware or software trigger. This provides a simple means of synchronizing the data collection between multiple MPS scanners when an IEEE-1588 architecture is not in place. It also provides a means of synchronizing MPS scanners with other equipment that may not be IEEE-1588 compatible. The MPS4264 has several options for external triggering. The external trigger can be used in one of three ways:

1. Hardware Frame Trigger
2. Hardware Scan Trigger
3. Software Frame Trigger

Options 1 & 3, the “frame trigger” options allow an external trigger (either hardware or software) to release a single frame of data. This means for a single trigger pulse (or command) the MPS will return a single frame of averaged data. Option 2, the “Hardware Scan Trigger” option means that if the trigger input level is “high” (5Vdc < voltage < 15Vdc) the MPS will enter scan mode where it will scan and continuously collect data based on its internal clock. The MPS will continue to scan at this rate until any of the following conditions are met:

- The trigger input level is switched to “low” (<5Vdc)
- FPS term is met
- The scan is manually stopped with the STOP command

The external hardware trigger voltage is input through pins 4 & 5 of the power connector. The trigger is 5-15Vdc, 6.5mA. A wiring diagram for the power/serial/trigger cable can be found in “Serial Communications” on page 14.

The external software trigger is a command sent to the MPS. The trigger can be one of two commands, either a TRIG command or a <tab>. More information on the TRIG variable can be found in “Set Trigger (TRIG)” on page 55. Information on the software trigger can be found in “Software Frame Trigger” on page 50.

SETTING UP AN EXTERNAL FRAME TRIGGER

To configure the MPS4264 to scan with an external frame trigger:

1. Set the variable “TRIG” to 1
2. Set the “RATE” variable to a rate that is approximately 25% faster than the intended trigger rate
3. Set the MPS into “scan” mode by sending the command: SCAN
4. Supply a trigger, either hardware trigger pulse or a

software trigger. For each trigger pulse, the MPS will return a single frame of data

SETTING UP AN EXTERNAL SCAN TRIGGER

To configure the MPS4264 to scan with an external scan trigger:

1. Set the variable “TRIG” to 2
2. Set the “RATE” variable to the desired data collection rate
3. Supply a “high” voltage (5Vdc<voltage<15Vdc). The MPS will collect data at the defined rate until the trigger input voltage is lowered below the 5Vdc trigger level.

AUTO SCAN OPTION

When TRIG is set to 3, the Auto Scan feature is enabled. This commands the MPS4000 to start scanning after the power up sequence has completed. The type of scan data output to the FTP server is set according to the format settings (CSV, ASCII or Binary). Setting TRIG to a 3 overrides the ENFTP setting because the only output available when auto scan is enabled is the FTP server.

When the FTP server is not ready to handle a connection from the MPS4264, the MPS will continue to attempt to connect until it succeeds, a STOP command is received from the serial interface or a STOP is received from the network interface. Stopping the auto scan by connecting with Telnet and issuing a STOP command is the preferred method used to change the TRIG setting from 3 when other trigger options are desired.

Once scanning is started, it continues until a STOP is received, as described above. Simply shutting off the MPS will stop the sending of data to the FTP server. Normally this will cause the FTP server to begin a time out period waiting for data. After the time out period, the FTP server should close the connection, flush any unsaved data to disk and close the file. The time out period of the FTP server is usually configurable. FTP servers may vary in their operation so check your FTP server manual for settings.

STATISTICAL SCAN OPTION

The MPS4264 includes an integrated function for analyzing the pressure data and providing statistical information on the data in real time. The MPS can output seven statistical parameters. They are: rolling average, rolling maximum, rolling minimum, rolling RMS, rolling standard deviation, rolling average excluding outliers at 3 sigma and number of overloads. The original data is also output. Because of the significant increase in data and processor load, the Statistical Scan option limits the overall data

output rate to 25Hz/channel with a buffer depth of 256.

All calculations exclude overloads where an overload is 150% of highest or lowest value in table. A list of the current overload value is found for each channel through the MXMN command (see “Max-Min Calibration Table Values” on page 49 for more information on the MXMN command.)

To configure statistical scan data a SET FORMAT of parameter B to S is required for data to be sent out the binary server port or parameter F to S is required for data to be sent out to the FTP server (see “Set Scanning Data Format (FORMAT)” on page 55 for more information on the FORMAT variable.) Statistical data is delivered to the web browser by clicking on the Statistical Scan link.

The depth of the subset FIFO is set via the SET OPTIONS <fast scan> <read mode> <sub set depth> variable in the LIST S group (see “Set Options (OPTIONS)” on page 56 for more information on the OPTIONS variable.) The valid range is from 2 to 256.

The packet format for a statistical packet is shown in “Table 5-6: Statistical Binary Data Format” on page 38.

When a client is connected to the binary server port, statistical scan data is sent out from the binary server port. A web browser must not be connected and the FORMAT variable must be set to S. The packet format is the same as shown in “Table 5-6: Statistical Binary Data Format” on page 38. Scanning may be started via the Telnet connection with either a SCAN or MSCAN command when a binary client is connected. Data is sent out the binary server port.

FAST SCAN OPTION (850Hz - 2500Hz)

Fast scan allows the MPS4264 to scan a reduced number of channels at a faster rate than is possible with all 64 channels. The set of reduced channels for a fast scan is 16 channels.

Not just any channel can be placed in the 16 channel fast scan group. Fast scan mode reduces the channel count down to one channel per A/D converter. Because of the association between channel numbers and A/Ds, the 16 channel fast scan group must contain 16 channels and start with channel 1, 2, 3 or 4. The channels that are part of each 16 channel fast scan group are shown in the table below.

To enable “Fast Scan” mode, the first term of the OPTIONS variable is set to the desired starting channel number. The setting of the first term of the OPTIONS variable will

correspond to the channel list in Table 3-4 below. The fast scan setting of 1, 2, 3 or 4 applies to the binary (B) and CSV (C) format settings. An error is produced if used with non B and C settings. An error is also produced if the “starting channel” is not 1, 2, 3 or 4. Format C is only useful for checkout and setup because the fast rate can not be maintained with Format C. Format B is the only setting that can maintain fast rates. The channels that are not part of the selected fast scan group have invalid data. External frame trigger can be used for fast scan. Starting a fast scan from the Web Browser is not supported.

More information on the OPTIONS variable can be found in “Set Options (OPTIONS)” on page 56. The data packet for fast scan data is defined in “Table 5-7: Fast Scan Binary Data Format” on page 39.

TABLE 3-4: FAST SCAN CHANNEL GROUPS

STARTING CHANNEL	NUMBER OF CHANNELS	CHANNELS IN SCAN
1	16	1,5,9,13,17,21,25,29,36,40, 44,48,52,56,60,64
2	16	2,6,10,14,18,22,26,30,35,3 9,43,47,51,55,59,63
3	16	3,7,11,15,19,23,27,31,34,3 8,42,46,50,54,58,62
4	16	4,8,12,16,20,24,28,32,33,3 7,41,45,49,53,57,61
0	64	All channels, 1-64

With “Fast Scan” enabled, the RATE variable is used to determine the data output rate of the limited group of channels. Data can be scanned on 1/4 of the channels at up to 2500Hz.

[INTENTIONALLY LEFT BLANK]

SECTION 4: HARDWARE

PNEUMATIC CONNECTIONS

The MPS modules has several different pneumatic connections. All pneumatic inputs are located on the top of the module. There are two basic groups of pneumatic connection types on the module: pressure measurement connection and configuration/support connections. All pneumatic inputs are through one of three removable headers for easy plumbing and system re-configuration. The measurement ports are located on the outer edge of the modules and are identified as 1-32 on one pneumatic header and 33-64 on the opposite pneumatic header. These 64 input ports are directly connected to one of the individual pressure transducers within the module. One the same pneumatic headers as the measurement input ports are the purge input ports. Each purge input port is used as the purge supply feed for purging the input lines on the respective side of the module.

The remaining configuration/support type pressure inputs are all located on the center removable header. They are: CAL CTL, REF, CAL and Px CTL. These ports are used for applying known pressures during calibration, cycling the internal valve between states and connecting the reference manifold to a known, stable ambient pressure.

The 64 measurement input ports (or Px ports) are available with 0.042" or 0.031" OD tubes. The remaining tubes for the PRG, CAL, REF and CTL inputs will always be 0.063" OD.

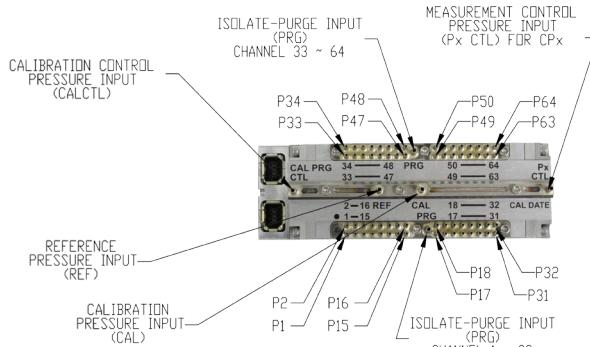


FIGURE 4-1: PNEUMATIC INPUTS, CPx VARIANT

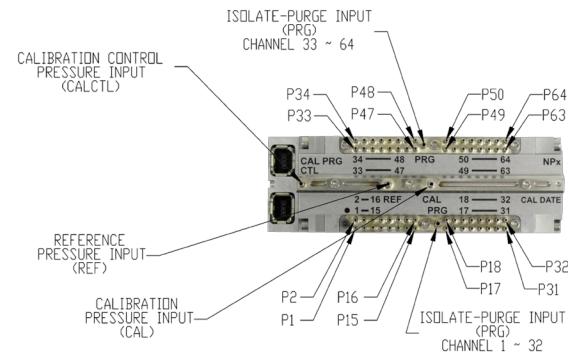


FIGURE 4-2: PNEUMATIC INPUTS, NPx VARIANT

BASICS OF CONNECTING INPUT TUBES

All pneumatic connections are 0.063", 0.042" or 0.031" OD stainless steel bulged tubulations. These tubes are designed to have a tight-fitting plastic tube slid over the OD of the tube. For all low pressure connections (50psi and less) the plastic tubing can simply be slid over the tubulation and the connection is complete. For higher pressure applications, when using Urethane or Vinyl tubing, a helical spring clamp over the OD of the plastic tubing is recommended. The clamp is slid over the tubing and located around the apex of the bulged of the tubulation. This helps hold the plastic tubing in place and prevent leaks.

Installing the plastic tubing over the stainless steel tubing can be done much easier by using Scanivalve's special "Tubing push-on tool" P/N: TPOTL-XXX. This tool is offered in a variety of sizes to work for tubing from 0.031" to 0.125" OD.

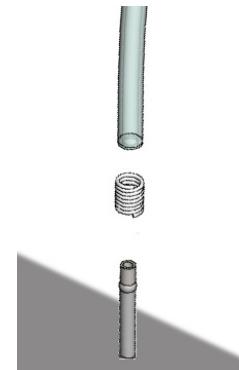


FIGURE 4-3: TUBULATION CONNECTION WITH HELICAL SPRING CLAMP

PRESSURE MEASUREMENT (Px) PORTS

Each MPS4264 scanner module has 64 pressure measurement ports, or “Px” ports. They are labeled as 1-64. Each of these ports are connected to a discrete pressure transducer. If any Px ports are not being used, it is recommended that they be plugged to prevent dust or any debris from clogging the port or contaminating the internal calibration valve.

CALIBRATION (CAL) PORT

The MPS module has a single “CAL” port in the center of the middle input header. The “CAL” port provides a means of applying a known calibration pressure to a single port and have it manifolded to the positive side of all of the transducers. During normal operations when the module is in measurement mode this port is internally blocked off. Any pressure applied to the “CAL” port will not reach the transducers unless the internal calibration valve has been configured to direct the calibration pressure to the transducers. For most low pressure applications (below 15psi) where accuracy is important, the “CAL” port should be connected to a known, stable static location when not being used to apply calibration pressures. This ensures that when a pneumatic zero offset calibration (CALZ) is performed no unwanted offsets are introduced.

REFERENCE (REF) PORT

The ‘REF’ port ties into a manifold that connects the back side (or negative side) or all transducers together. During most applications, the ‘REF’ port of low pressure modules (below 15 psi) should be routed to a known, stable static location. Typically this “reference” location will be a wind tunnel static port, a static barrel or in flight test applications the aircraft static system. This ensures that when a zero offset calibration (CALZ) is performed no unwanted offsets are introduced. During a calibration, positive pressures will be applied through the ‘REF’ port to perform the negative portion of the calibration.

For applications where an elevated reference pressure will be used, the maximum pressure applied to the REF input is 50psi.

PURGE (PRG) PORTS

The purge supply (PRG) port is the input for the pressure that will be used to purge the pressure measurement (Px) lines. ***The supplied purge pressure can be up to 75psi, and must always be clean, dry instrument grade air or nitrogen.*** In order for the purge supply pressure to successfully feed to and clear the Px input lines, the MPS must be in Calibration/Purge mode. This is done by applying control pressure to the “CAL-PRG CTL” port.

Recommended purge pressures are based off of the module’s pressure range and are as follows:

Module Pressure Range	Recommend Purge Pressure
4" H ₂ O	1.5 PSI
8" H ₂ O	3.0 PSI
1 PSI	10 PSI
5 PSI	25 PSI
15 PSI	75 PSI

CALIBRATE MODE CONTROL (CAL-PRG CTL) PORT

The “CAL-PRG” control pressure input is used to shift the internal calibration valves into Calibration/Purge mode. In this mode, the positive side of all transducers are tied together and connected to the CAL input, and the individual Px input tubes are manifolded together internally and connected to the purge supply (“PRG”) ports. ***The maximum input pressure for the CAL-PRG CTL input is 120psi.*** Clean, dry instrument grade air or nitrogen should be used.

For modules that are configured with the “NPx” actuation, 90-120psi is recommended to cycle the valve completely and reliably.

For modules that are configured with the “CPx” actuation option, the required control pressure to cycle the valve is 60-70psi. Often this can be done with less pressure, but 60psi is recommended for reliable operation of the valve. For “CPx” modules, in normal, low-vibration environments the control pressure does not need to be continuously supplied to this port to hold the valve in the Calibrate mode. Applying the pressure for 2-4 seconds will cycle the valve into the calibrate mode where it will stay until cycled out using the Px CTL port. For high vibration environments, or where it is absolutely critical that the module be held in calibrate mode, it is recommended to continuously supply the pressure to this port to hold the valve in the desired state.

MEASUREMENT MODE CONTROL (Px CTL) PORT

The “Px CTL” control pressure input is used to shift the internal calibration valve into the “Measurement” mode. In this mode, each of the individual transducers are tied to the discrete input (Px) tubes. The CAL and PRG input ports are blocked off. This input port is only present on modules configured with the “CPx” actuation option. ***The maximum input pressure for the Px CTL input is 70psi.*** To actuate or cycle the valve 60-70psi is recommended. Often this can be done with less pressure, but 60psi is recommended for reliable operation of the valve. Clean, dry instrument grade air or nitrogen should be used.

In normal, low-vibration environments the control pressure does not need to be continuously supplied to this port to hold the valve in the measurement mode. Applying the pressure for 2-4 seconds will cycle the valve into the measurement mode where it will stay until cycled out using the Px CTL port. For high vibration environments, or where it is absolutely critical that the module be held in measurement mode, it is recommended to continuously supply the pressure to this port to hold the valve in the desired state.

REMOVABLE PRESSURE INPUT HEADERS

All pneumatic inputs to the MPS are through three removable input headers. Each of the three headers is held in place with three #1-72 screws. A 1/16" Allen/hex wrench fits the screws.

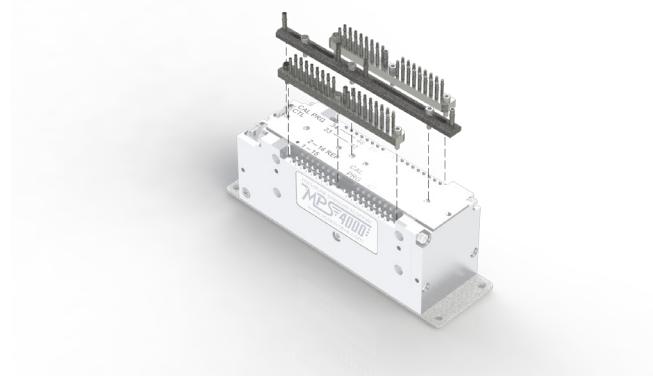


FIGURE 4-4: REMOVABLE INPUT HEADERS

The center header is polarized and cannot be installed backwards. However, the two Px input headers are interchangeable and can be installed on the opposite sides. To keep the manifolds unique and identified, the header for inputs 1 through 32 has a single black oxide screw, whereas all of the rest are nickle plated and silver in color. The "correct" location for the black oxide screw is identified on the MPS with a solid black dot. See Figure 4-5 below:

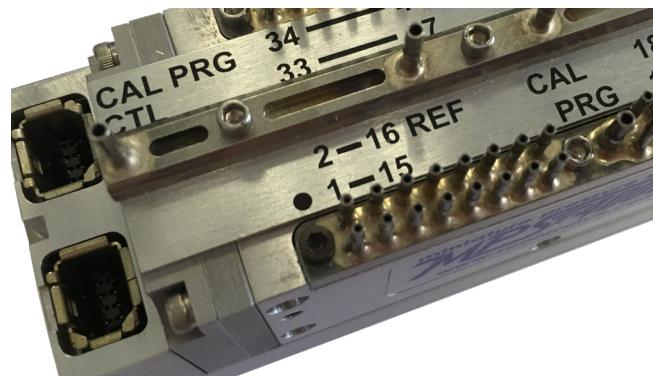


FIGURE 4-5: HEADER IDENTIFIER MARK

Be aware, that when the center header is removed the

identification plate will no longer be held in place. Like the center header it is indexed and cannot be installed backwards, but should be handled with care. The O-rings supported in the identification plate should be kept clean and watched to make sure they stay in place. Replacement or additional headers can be purchase separately as needed.

CALIBRATION VALVE

The MPS4264 module has two mirrored calibration valves, each servicing 32 inputs. It is a slider-style valve with the aluminum shuttle riding on ball bearings. The O-rings that move with the valve shuttle ride on a proprietary plastic "bearing plate" that minimizes friction to reduce the actuation force and the wear on the O-rings.

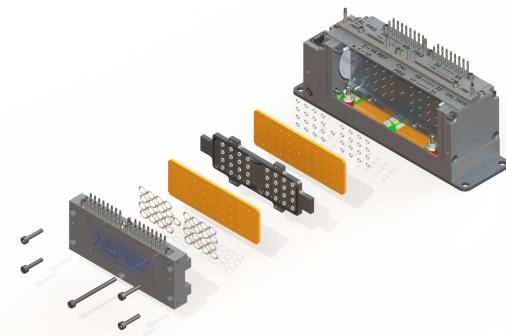


FIGURE 4-6: VALVE EXPLODED VIEW

The calibration valve has two physical states, "Measurement" and "Calibrate/Purge." In measurement mode, measurement pressures from the individual input tubes is allowed to pass directly through the shuttle valve to the pressure sensors. In calibrate/purge mode, the measurement pressures from the individual input tubes are blocked off and isolated from the pressure sensors. Instead, the "CAL" input is manifolded and connected to the positive side of all transducers. Additionally, the two "Purge" input lines are connected to the input tubes on their respective valves allowing the user to clear debris or moisture from the input lines by blowing from the module back out to the test article.



CAUTION! Do not cycle the valve with pressure supplied to either the CAL or PRG inputs.

It is uncommon, but possible, that while the valve is

switching states (between measurement and calibration/purge modes) that any pressure applied to the "CAL" or the "PURGE" ports could "sneeze" past the sealing O-rings and reach the pressure sensor. To prevent risking potential damage to the sensors, the valve should always be cycled with no "CAL" or "PRG" pressure supplied.

Also note that when the valve is in the calibrate/purge state, all 32 measurement inputs on each valve are tied together and connected to the "PRG" input. For most applications this is not a problem, but it will allow "mixing" of measurement pressures between channels. In cases where humid or dirty air is being measured, this "mixing" would provide a flow path in a higher pressure input and back out of a lower pressure input, carrying moisture or debris with it. If it is the case where the module will be help in calibrate/purge mode for an extended period of time while still connected to dirty or humid measurement pressures, a purge pressure great enough to prevent flow from the test article to the scanner should be supplied. For most applications however, this is not needed.

ELECTRICAL CONNECTORS

The MPS4264 has two electrical connectors located adjacent to each other on one end of the module. Both connectors are TE Connectivity "Mini I/O" series and are polarized to prevent improper connections. The connectors are latching to ensure a reliable connection. To disconnect the connectors, press down on the latch (shown in Figure 4-7 below) and pull up on the connector.

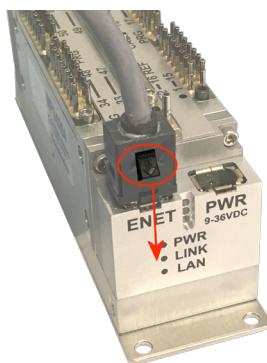


FIGURE 4-7: POWER/ETHERNET CONNECTOR LATCH

POWER CONNECTOR & POWER REQUIREMENTS

The power connector servers three purposes. Primarily it is the power connector, bringing in the 9-36Vdc required to power the MPS. Additionally, it serves as the external trigger input (pins 4 & 5) as well as the serial communications connection (pins 1-3). See "Appendix F - Cable Drawings and Pinouts" on page 88 for cable drawings and pinout.

The MPS accepts 9-36Vdc and has a maximum power consumption of 3.5W.



CAUTION! Do not make or break the power connector with power applied! Doing so risks damage to the module!

ETHERNET CONNECTOR


The Ethernet connection is the primary means of communications with the MPS4264. The connection does support MDIX Auto-Crossing. See "Appendix F - Cable Drawings and Pinouts" on page 88 for cable drawings and pinout.

MPSTCU-THERMAL CONTROL UNIT

The MPS (Miniature Pressure Scanner) line of Thermal Control Units is available for applications where temperature changes may be great enough to exceed the temperature compensated range of the scanner. Exceeding the compensated temperature range can induce errors in the pressure measurements. See "Section 1: Specifications" on page 7, for more information on the compensated temperature range. MPSTCU applications include flight test, automotive, wind turbine, wind tunnel, and engine tests. When used in a flight test environment, it is recommended that the TCU case be vented.

Each MPSTCU includes a rugged IP-54 rated aluminum enclosure, pneumatic connectors, and a power/data connection. The MPSTCU offers an optional heater circuit which can provide a controlled temperature environment while the MPSTCU is subject to a cold environment (up to -60°C). Additional, the MPSTCU also offers an optional cooling port to provide cooling air when the MPSTCU is subject to temperatures from 60-120°C. When the Cooling Kit is used, 3.0CFM of 23°C air is required while operating the MPSTCU in a 125°C environment.

For wiring diagrams of MPSTCU cabling, see "MPS4264TCU Cable" on page 90.

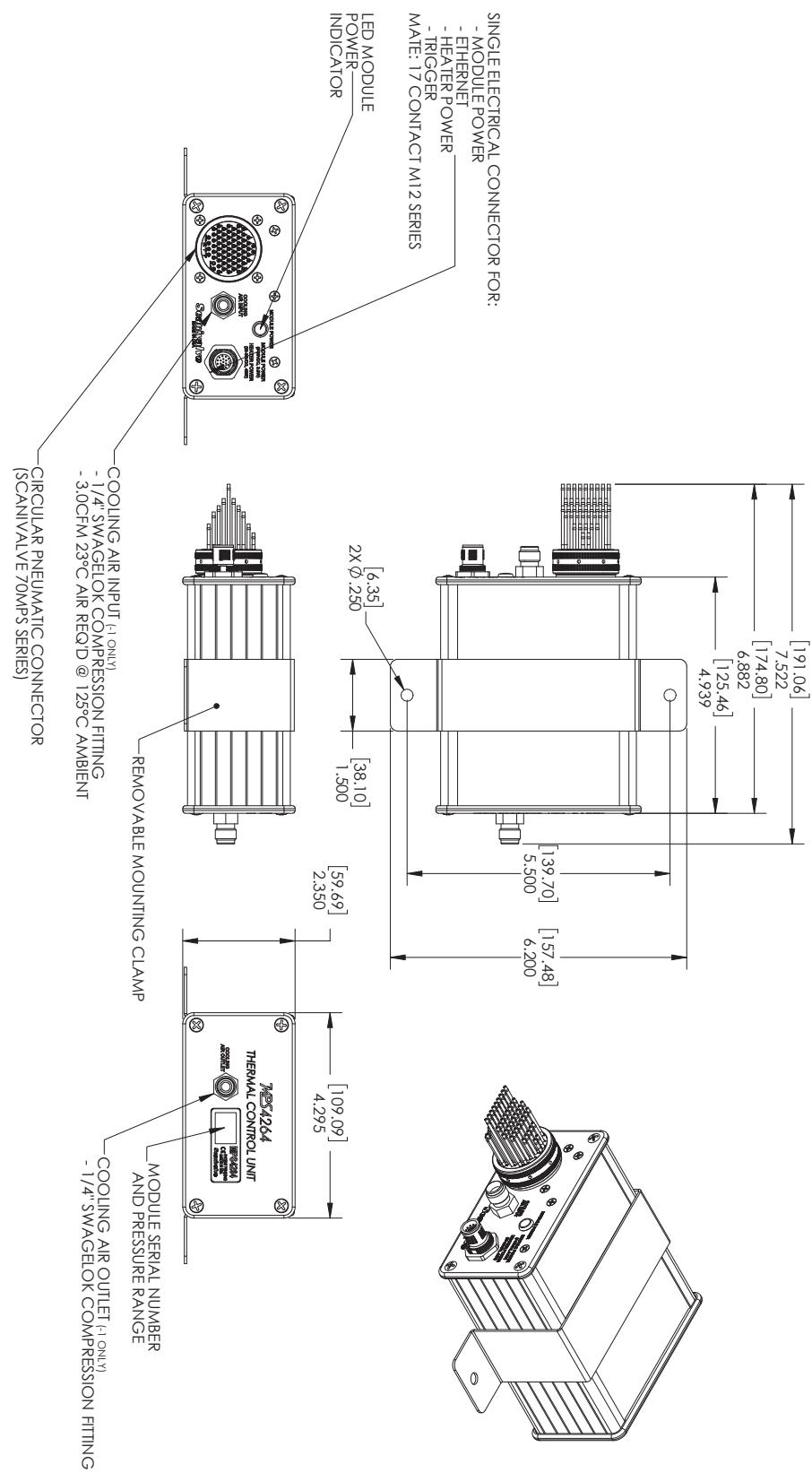MPS4264TCU OVERVIEW DRAWING

[INTENTIONALLY LEFT BLANK]

SECTION 5: SOFTWARE

OVERVIEW

The MPS4264 is essentially a single user network server that supports connections from a Telnet Client, Web Browser Client, Binary Client, FTP Client (for data transfers) and FTP Server (for file transfers). The MPS4264 supports serial communication that is solely used to set a limited group of configuration variables for diagnostics or when the IP address is not known. The MPS4264 also supports a limited number of commands from a multicast master. It also acts as a multicast master. This design allows a single MPS4264 to reissue commands to multiple MPS4264s acting as slaves.

COMMAND FORMAT

Each of the commands are explained with the following sections: command, syntax, arguments, description, and returns.

COMMAND lists the name of the command.

SYNTAX lists the format of the command. The following conventions are used:

BP - Boldface letters indicate command keywords and operators. Within the discussion of syntax, bold type indicates that the text must be entered exactly as shown.

<parameter> - Words in <> indicate place holders for parameters you must supply.

[parameter] - Words in [] indicate place holders for optional parameters that you are not required to supply.

data - Example data is shown in *italics*.

data - Example commands are show in **bold**.

, - Commas separate options, only one of the options may be used.

<CR> - The carriage-return key, sometimes marked as a bent arrow, Enter, or Return on the key board, is called <CR>.

Spaces, as used in the syntax, are entered as spaces.

DESCRIPTION describes the function of the command.

RETURNS lists the format of the information that the unit returns to the host.

A PROMPT (>) will be output when the MPS4264 is ready to accept a command.

TCP/IP does not guarantee that packet boundaries will be maintained between a Host and a MPS4264. Therefore, ALL commands from a Host MUST be terminated properly with one of two options using the NL configuration variable. The

two options are:

CR-LF (ASCII 13 - ASCII 10) or CR (ASCII 13)

The maximum string length for commands send to the MPS4264 is 79 characters. Any command received that is longer than 79 characters will be discarded and an error will be generated.

When a communications variable is modified (those in the LIST IP group), the MPS4264 program must be restarted, preferably with the RESTART command, in order for the changes to take effect.

OUTPUT SCAN DATA FORMATTING AND ROUTING

The MPS4000 supports four data destinations. They are: FTP Server, Telnet Client, Web Browser Client, and Binary Client. The data destination is based on the connections, configurations, and source of command. The following table describes the output data routing based on the factors defined in "Table 5-1: Data Destinations" on page 34.

TABLE 5-1: DATA DESTINATIONS

CONFIGURATION	DESTINATION NATION	SUPPORTED FORMATS									
		WB	Bin	TN	FTP	UDP	FTP & UDP	FTP	UDP	FTP & UDP	None
Web Browser Connected	•										
Binary Server Connected (SVRSEL 2)											
Telnet Connected											
FTP Enabled											
UDP Enabled (SVRSEL 3)											
Auto Scan Enabled											
Scan Command From Telnet											
Scan Command From Multicast Master (SVRSEL 1)											
Scan Command From Serial											
Scan Command From External Trigger											
Data To:	•	•	•	•	•	•	•	•	•	•	•
ASCII											
CSV											
VT100											
Binary LabVIEW											
UDP Binary Packet											
Binary Packet											
Binary Statistical Packet											
Web Browser Bar Graph											

Notes:

Any new connection that occurs while scanning will cause scanning to stop.

Blank cells indicate "Not Possible"

Bin Con and Multi Server are multiplexed and cannot both be active at the same time.

When the web browser is connected all commands are ignored from serial and Telnet with an error message

WB – Web Browser, Bin – Binary Server, TN – Telnet, FTP – File Transfer Protocol Server

ASCII DATA FORMATS

Scan data from the MPS4000 can be formatted in six different ways: ASCII (A), Comma Separated Variable (CSV or C), formatted for VT100 terminals (F), Binary (B) LabVIEW Binary (L) and Statistical Binary (S). The format is selected by setting the FORMAT configuration variable to A, C, F, B, L, or S respectively. It can also be selected by clicking on the appropriate link on the “Scan” page when using the Web Browser.

The FORMAT configuration variable sets the format for Telnet (T), FTP (F) and Binary Server (B) separately. See the description of FORMAT for detailed description.

Data output in ASCII, CSV and FORMATTED are typically used for low speed data transfers, less than 10 Hz. Binary data output is used for high speed data transfer up to 850 Hz. Appendix A shows the binary packet format and Appendix B shows the statistical binary packet format.

BASIC ASCII DATA OUTPUT

“Table 5-2: ASCII Data Output” shows an example of ASCII data output (SET FORMAT T A):

TABLE 5-2: ASCII DATA OUTPUT

1 1 0.0021 30.44
1 2 -0.0009 29.88
1 3 -0.0010 30.31
1 4 -0.0015 30.13
1 5 0.0020 30.31
1 6 -0.0004 30.06
1 7 0.0006 29.94
1 8 -0.0012 30.19
1 9 0.0023
1 10 0.0006
1 11 -0.0009
.....
1 59 -0.0004
1 60 0.0021
1 61 -0.0005
1 62 -0.0007
1 63 -0.0003
1 64 0.0021

Notes:

Configuration selected with by setting the “FORMAT” variable to “A”.

Column 1 is the frame number

Column 2 is the channel number

Column 3 is the pressure in either engineering units or raw counts.

Column 4 for the first 8 rows is the temperature of the 8 temperature sensors in degrees C

FORMATTED ASCII DATA OUTPUT

"Table 5-3: Formatted ASCII Data Output" shows an example of ASCII data output formatted for a VT100 Terminal (SET FORMAT T F):

TABLE 5-3: FORMATTED ASCII DATA OUTPUT

Frame= 1
T1= 35.1 C T2= 34.7 C T3= 35.1 C T4= 34.9 C
T5= 35.1 C T6= 34.8 C T7= 34.7 C T8= 34.9 C
01= 0.0028 02= -0.0016 03= -0.0011 04= -0.0014 05= 0.0020 06= -0.0003
07= -0.0003 08= -0.0014 09= 0.0017 10= 0.0011 11= -0.0002 12= -0.0009
13= 0.0032 14= 0.0015 15= 0.0028 16= -0.0005 17= 0.0018 18= 0.0022
19= -0.0006 20= -0.0019 21= -0.0000 22= -0.0007 23= -0.0018 24= -0.0022
25= 0.0029 26= 0.0001 27= 0.0000 28= -0.0006 29= -0.0000 30= 0.0010
31= 0.0004 32= -0.0007 33= -0.0005 34= -0.0018 35= 0.0008 36= 0.0011
37= 0.0003 38= -0.0001 39= -0.0001 40= 0.0029 41= -0.0005 42= -0.0016
43= -0.0003 44= 0.0019 45= 0.0003 46= -0.0002 47= 0.0001 48= 0.0021
49= -0.0002 50= 0.0003 51= 0.0004 52= 0.0016 53= 0.0014 54= 0.0002
55= 0.0009 56= 0.0007 57= -0.0006 58= -0.0010 59= -0.0009 60= 0.0024
61= -0.0008 62= -0.0006 63= 0.0002 64= 0.0026
Notes:
Configuration selected with by setting the "FORMAT" variable to "F".
This format should only be selected for displaying data on a screen.
The 8 temperatures are shown as T1 through T8.
The 64 pressures are shown as 01 through 64
This format uses VT-100 control codes to set cursor positions for column alignment.

ASCII CSV DATA OUTPUT

"Table 5-4: CSV Data Output" shows an example of Comma Separated Variable data formatted (SET FORMAT T C):

TABLE 5-4: CSV DATA OUTPUT

Frame,Seconds,Tx1,Tx2,Tx3,Tx4,Tx5,Tx6,Tx7,Tx8,Px1,Px2,Px3,Px64
1,0.100,40.38,39.88,40.38,40.13,40.31,40.00,39.88,40.13,0.0030,-0.0011,-0.0010,-0.0025
2,0.200,41.50,41.00,41.50,41.25,41.44,41.13,41.00,41.25,0.0022,-0.0013,-0.0014,-0.0024
3,0.300,41.50,41.00,41.50,41.19,41.44,41.13,41.00,41.25,0.0017,-0.0017,-0.0006,-0.0030
4,0.400,41.50,41.00,41.50,41.19,41.44,41.13,41.00,41.25,0.0013,-0.0018,-0.0005,-0.0026
5,0.500,41.50,41.00,41.50,41.19,41.44,41.13,41.00,41.25,0.0021,-0.0013,-0.0005,-0.0028
6,0.600,41.50,41.00,41.50,41.19,41.44,41.13,41.06,41.25,0.0025,-0.0015,-0.0009,-0.0031
7,0.700,41.50,41.00,41.50,41.19,41.44,41.13,41.06,41.25,0.0027,-0.0011,-0.0015,-0.0026
8,0.800,41.56,41.00,41.50,41.19,41.44,41.13,41.06,41.25,0.0026,-0.0004,-0.0015,-0.0025
Notes:
Configuration selected with by setting the "FORMAT" variable to "C".
The output shown omits channels 4 through 63 for simplicity
Each column is separated by a comma and each row is delimited by a CR LF
This format is easily imported to Excel or Open Office spread sheets.

BINARY DATA FORMATS

Binary data formats follow network byte order (big endian).

BINARY DATA FORMAT**TABLE 5-5: BINARY DATA FORMAT**

FUNCTION	BYTES	DATA TYPE	DESCRIPTION
Packet Type	4	Integer	A Hex
Packet Size	4	Integer	Size in bytes of this packet (348)
Frame Number	4	Integer	The current frame number if in the scan mode
Scan Type	4	Integer	0-Neg, 1-Pos, 2-A/C
Frame Rate	4	Float	Scanning rate in Hz
Valve Status	4	Integer	0-Px, 1-Cal
Units index	4	Integer	Units conversion index (add table)
Units conversion factor	4	Float	Conversion factor from PSI to selected units
PTP Scan Start Time (sec)	4	Integer	Scan start time in seconds
PTP Scan Start Time (ns)	4	Integer	Scan start time in nanoseconds
External Trigger Time	4	U Integer	External Trigger Time in us
Temperatures (8 temperatures)	32	Floats	Array of 8 temperature values as read from the 8 temperature chips.
Pressures	256	Float or Integer	Array of 64 pressure values. Float if EU is selected, Integer is RAW is selected. The float or integer is selected based on the units.
Frame Time	4	Integer	Time the frame occurred (seconds)
Frame Time	4	Integer	Time the frame occurred (nanoseconds)
External Trigger Time	4	Integer	Time the external trigger occurred (seconds)
External Trigger Time	4	Integer	Time the external trigger occurred (nanoseconds)

STATISTICAL BINARY DATA FORMAT

TABLE 5-6: STATISTICAL BINARY DATA FORMAT

FUNCTION	BYTES	DATA TYPE	DESCRIPTION
Packet Type	4	Integer	11Hex
Packet Size	4	Integer	Size in bytes of this packet (2140)
Frame Number	4	Integer	The current frame number if in the scan mode
Scan Type	4	Integer	0-Neg, 1-Pos, 2-A/C
Frame Rate	4	Float	Scanning rate in Hz
Valve Status	4	Integer	0-Px, 1-Cal
Units index	4	Integer	Units conversion index (add table)
Units conversion factor	4	Float	Conversion factor from PSI to selected units
PTP Scan Start Time (sec)	4	Integer	Scan start time in seconds
PTP Scan Start Time (ns)	4	Integer	Scan start time in nanoseconds
External Trigger Time	4	U Integer	External Trigger Time in us
Temperatures (8 temperatures)	32	Floats	Array of 8 temperature values as read from the 8 temperature chips.
Pressures	256	Float or Integer	Array of 64 pressure values. Float if EU is selected, Integer is RAW is selected. The float or integer is selected base on the units.
Frame Time	4	Integer	Time the frame occurred (seconds)
Frame Time	4	Integer	Time the frame occurred (nanoseconds)
External Trigger Time	4	Integer	Time the external trigger occurred (seconds)
External Trigger Time	4	Integer	Time the external trigger occurred (nanoseconds)
Rolling Average	256	Float	Array of rolling average of pressures
Max Value	256	Float	Array of rolling maximum pressure
Min Value	256	Float	Array of rolling minimum pressure
RMS Value	256	Float	Array of rolling RMS of pressure
Standard Deviation	256	Float	Array of rolling standard deviation of pressure
Rolling Average Excluding Outliers	256	Float	Array of rolling average of pressure excluding outliers at 3 sigma

FAST SCAN BINARY DATA FORMAT**TABLE 5-7: FAST SCAN BINARY DATA FORMAT**

FUNCTION	BYTES	DATA TYPE	DESCRIPTION
Packet Type	4	Integer	10Hex
Packet Size	4	Integer	Size in bytes of this packet (348)
Frame Number	4	Integer	The current frame number if in the scan mode
Scan Type	4	Integer	0-Neg, 1-Pos, 2-A/C
Frame Rate	4	Float	Scanning rate in Hz
Valve Status	4	Integer	0-Px, 1-Cal
Units index	4	Integer	Units conversion index (add table)
Units conversion factor	4	Float	Conversion factor from PSI to selected units
PTP Scan Start Time (sec)	4	Integer	Scan start time in seconds
PTP Scan Start Time (ns)	4	Integer	Scan start time in nanoseconds
External Trigger Time	4	U Integer	External Trigger Time in us
Temperatures (8 temperatures)	32	Floats	Array of 8 temperature values as read from the 8 temperature chips.
Pressures	256	Float or Integer	Array of 64 pressure values. Float if EU is selected, Integer is RAW is selected. The float or integer is selected based on the units. When fast scan is enabled only 16 of the 64 channels are valid.
Frame Time	4	Integer	Time the frame occurred (seconds)
Frame Time	4	Integer	Time the frame occurred (nanoseconds)
External Trigger Time	4	Integer	Time the external trigger occurred (seconds)
External Trigger Time	4	Integer	Time the external trigger occurred (nanoseconds)

FILE NAMES AND FORMATS

The MPS4264 can be configured to output a variety of different files types and configurations. “Table 5-8: File Name Extensions” below shows the file name extensions used based on the setting of the FORMAT variable.

TABLE 5-8: FILE NAME EXTENSIONS

SETTING OF FORMAT VARIABLE	TYPE OF SCAN	FILE EXTENSION
SET FORMAT F A	ASCII	nnnn.txt
SET FORMAT F B	Binary	nnnn.dat
SET FORMAT F C	Comma Separated Variable	nnnn.csv
SET FORMAT F S	Statistical Scan	nnnn.dat

GENERAL CONTROL COMMANDS**BOOT LOADER VERSION**

DESCRIPTION	Shows the current software version of the boot loader.
SYNTAX	BLVER
RETURNS	The boot loader version followed by the prompt when ready for the next command.
EXAMPLE	BLVER <i>BOOTLOADER VERSION: Ver 2.03</i> >

FIRMWARE FILES CHECK

DESCRIPTION	Checks the current files loaded on the MPS against the index in the .TAR file. This verifies file names and creation dates.
SYNTAX	TARCHK
RETURNS	List of files and error summary
EXAMPLE	TARCHK <i>Starting TARCHK</i> <i>Testing file MPS4000.BIT</i> <i>Testing file Mps4000.hex</i> <i>Testing file calval.htm</i> <i>Testing file databar.htm</i> <i>Testing file favicon.ico</i> <i>Testing file home.gif</i> <i>Testing file home.jpg</i> <i>Testing file home.htm</i> <i>Testing file idframe.htm</i> <i>Testing file index.htm</i> <i>Testing file navframe.htm</i> <i>Testing file remote.htm</i> <i>Testing file setup.htm</i> <i>Testing file svlogo.gif</i> <i>Testing file term.htm</i> <i>Finished TARCHK with 0 errors</i>

FORMAT FLASH DISK

DESCRIPTION	Formats the flash memory of the MPS. WARNING: Formatting the flash memory is considered to be an advanced operation. Use caution before proceeding.
SYNTAX	FDISK
RETURNS	<pre>>fdisk Type FDISKCONFIRM to confirm FDISK or STOP to escape</pre> <pre>>fdiskconfirm Start format Erasing device 100 Done device 100</pre>
EXAMPLE	<pre>>fdisk Type FDISKCONFIRM to confirm FDISK or STOP to escape</pre> <pre>>fdiskconfirm Start format Erasing device 100 Done device 100</pre>

MULTIPLE UNIT FIND

DESCRIPTION	Lists all of the MPS4264 devices in the multicast cluster of the master.
SYNTAX	MFIND
RETURNS	A list of MPS modules found in the multicast “cluster” followed by a prompt ready for the next command.
EXAMPLE	MFIND <i>Found device 0 IP Address 191.30.90.9</i>

MULTIPLE UNIT STOP

DESCRIPTION	Sent to stop multiple units when working with multiple devices. When the master receives the MSTOP command it issues a STOP command to the multicast address set in the MCAST variable. Each MPS member of this multicast cluster will stop scanning.
SYNTAX	MSTOP or <ESC>
RETURNS	Prompt ready for the next command.
EXAMPLE	MSTOP >

REBOOT

DESCRIPTION	Restarts the MPS4000.
SYNTAX	REBOOT or RESTART
RETURNS	Nothing
EXAMPLE	REBOOT

STATUS

DESCRIPTION	Displays the current mode of the MPS4000
SYNTAX	STATUS
RETURNS	One of the following: READY – Ready mode SCAN – Currently scanning CAL – Performing a calibration VAL – Performing a validation CALZ – Performing a zero calibration CALM – Performing a manual calibration Followed by a prompt when ready for the next command.
EXAMPLE	STATUS <i>STATUS: READY</i> >

STOP

DESCRIPTION	Cancels all commands and returns the MPS4000 to Ready mode.
SYNTAX	STOP or <ESC>
RETURNS	Prompt ready for the next command.
EXAMPLE	STOP >

TAR FILE UNPACK

DESCRIPTION	Commands the MPS to unpack the designated .TAR file.
SYNTAX	TAR <filename>
RETURNS	List of files unpacked.
EXAMPLE	TAR MPS123.tar <i>Starting to de-tar files from fileMPS105.tar</i> <i>Extracting file MPS4000.BIT 340701 2016/05/18 18:29:30</i> <i>Extracting file Mps4000.hex 695026 2016/06/13 20:54:24</i> <i>Extracting file calval.htm 19639 2015/04/17 15:05:33</i> <i>Extracting file databar.htm 18158 2016/01/12 18:56:09</i> <i>Extracting file favicon.ico 1150 2014/03/14 15:19:55</i> <i>Extracting file home.gif 81414 1980/01/01 19:00:00</i> <i>Extracting file home.jpg 16814 2014/05/12 19:14:03</i> <i>Extracting file home.htm 761 2015/05/22 22:51:24</i> <i>Extracting file idframe.htm 3933 2015/06/01 17:08:13</i> <i>Extracting file index.htm 2067 2015/05/26 14:47:55</i> <i>Extracting file navframe.htm 2944 2015/05/27 18:36:18</i> <i>Extracting file remote.htm 8141 2015/04/16 19:48:57</i> <i>Extracting file setup.htm 10312 2016/03/08 16:40:52</i> <i>Extracting file svlogo.gif 2472 2015/06/01 17:03:21</i> <i>Extracting file term.htm 3101 2015/06/03 19:33:35</i> <i>Done extracting files from MPS105.tar</i> <i>Creating TAR.NDX file</i> <i>Done creating TAR.NDX file</i> >
NOTES	Capitalization of the filename entered is case sensitive and must match the file name exactly.

VALVE STATE

DESCRIPTION	Displays the current position of the Pneumatic Configuration Valve
SYNTAX	VALVESTATE
RETURNS	One of the following: CAL – Calibration mode PX – PX mode Followed by a prompt when ready for the next command.
EXAMPLE	VALVESTATE <i>PX</i> >

VERSION

DESCRIPTION	Shows the current software version
SYNTAX	VER
RETURNS	The version followed by the prompt when ready for the next command
EXAMPLE	VER <i>MPS Scanivalve (c) 2016 Ver 1.50</i> >

PTP (PRECISION TIME PROTOCOL) COMMANDS

CLEAR ACCUMULATED PTP STATISTICAL DATA

DESCRIPTION	Clears the accumulated data in the histogram buffers and restarts sampling.
SYNTAX	PTPCLEAR
RETURNS	A prompt when ready for the next command.
EXAMPLE	PTPCLEAR <i>PTP Histogram Cleared</i> >

DISPLAY PTP STATISTICS

DESCRIPTION	Displays the accumulated data in histogram form.
SYNTAX	PTPSTAT
RETURNS	OFM, INTERVAL, <interval value>, <interval -6>, <interval -5>, <interval -4>, <interval -3>, <interval -2>, <interval -1>, <interval 0>, <interval 1>, <interval 2>, <interval 3>, <interval 4>, <interval 5>, <interval 6> A prompt when ready for the next command.
EXAMPLE	PTPSTAT <i>OFM, INTERVAL, 100000, 114, 13, 10, 23, 21, 15, 36, 6, 12, 9, 23, 10, 129</i> >

GET PTP TIME

DESCRIPTION	Gets the current PTP time as used by the MPS. Time is adjusted by UTCOFFSET.
SYNTAX	GETTIME
RETURNS	Current Time yyyy/mm/dd hh:mm:ss sec nnnn ns
EXAMPLE	GETTIME <i>Current Time 2014/2/6 0:23:52.0 sec 3335 ns 439809720</i> >

GET SCAN START TIME

DESCRIPTION	Gets the start time of the last scan. This command is valid in all modes. If no scan has been run since power up the following time is returned: 2015/1/1 0:0:0.000000
SYNTAX	SST
RETURNS	yyyy/mm/dd hh:mm:ss.nnnnnnnn
EXAMPLE	SST <i>2016/4/29 8:45:25.000000</i> >

GET UTC OFFSET

DESCRIPTION	Gets the current UTC offset.
SYNTAX	GETUTCO
RETURNS	- The current difference between TAI time and UTC time, in seconds, as supplied by the local grandmaster. - UTC updated flag - Set when UTC offset is valid.
EXAMPLE	GETUTCO <i>Current UTC Offset 0 0 0</i> >

SET PTP TIME

DESCRIPTION	Sets the current PTP time as used by the MPS. Time is adjusted by UTCOFFSET. When a grandmaster is active on the network, the time set by the SETTIME command is overwritten by the grand master time.
SYNTAX	SETTIME yyyy/mm/dd hh:mm:ss
RETURNS	Current Time yyyy/mm/dd hh:mm:ss sec nnnn ns
EXAMPLE	SETTIME 2016/6/22 14:20:0.00 <i>Entered 2016/6/22 14:20:0</i> <i>Time is 1466605200</i>

SCAN AND CALIBRATION COMMANDS

Reference Scanivalve Document: MPS Field Val if performing a single temperature validation with an SPC4000.

CALIBRATION

DESCRIPTION	Starts a calibration for the MPS4264.
SYNTAX	<p>CAL [<F or FM>] Where: *F - Initiates a single temperature calibration, with no oven connection. Connects to the calibrator. *FM- Initiates a manual single unit single temperature calibration with no oven connection. Does not connect to the calibrator When no parameter is present, the unit performs a self calibration based on the calibration group settings. *Not yet released</p>
RETURNS	Calibration data
EXAMPLE	CAL >

MULTIPLE UNIT CALIBRATION

DESCRIPTION	When the master receives the MCAL command it starts the multiple device calibration sequence described in the section on multiple device calibration and validation.
SYNTAX	<p>MCAL [<F or FM>] Where: *F- Initiates a multiple unit single temperature calibration with no oven connection. Connects to the calibrator. *FM- Initiates a manual multiple unit single temperature calibration with no oven connection. Does not connect to the calibrator When no parameter is present, the unit performs an automated validation based on the calibration group settings *Not yet released</p>
RETURNS	Calibration data
EXAMPLE	MCAL >

CALIBRATION & VALIDATION

DESCRIPTION	Starts a calibration followed by a validation and places the MPS4264 into CAL mode followed by VAL mode. Reference Scanivalve Document: MPS Field Val if performing a single temperature validation with an SPC4000.
SYNTAX	CALVAL
RETURNS	Validation data
EXAMPLE	CALVAL >

MULTIPLE UNIT CALIBRATION & VALIDATION

DESCRIPTION	When the master receives the MVAL command it starts the multiple device calibration and validation sequence described in the section on multiple device calibration and validation. Reference Scanivalve Document: MPS Field Val if performing a single temperature validation with an SPC4000.
SYNTAX	MCALVAL
RETURNS	Validation data
EXAMPLE	MCALVAL >

CURRENT TEMPERATURE PLANE

DESCRIPTION	Returns the current temperature plane for the channel.
SYNTAX	CTP <channel> [<temp>] Where channel is a channel number from 1 to 64. Temp is the optional temperature. When omitted the current temperature is used.
RETURNS	Number of pressure points, channel number, temperature used, number of temperature used for average, temperature gain used to determine the pressure index, and temperature offset used to determine the pressure index. The channel, current temperature, pressure index, pressure, counts, M and B terms for the current or passed temperature.
EXAMPLE	CTP 1 10 CTP NumPress 21 Channel 1 Temp 10.00 0 Mndx 5.355829e-06 Bndx 14.198914 C 1 T 10.00 P 0 Prs -5.499897 Cnt -2708841 M 1.584424e-06 B -1.299782 C 1 T 10.00 P 1 Prs -4.950220 Cnt -2513036 M 1.663405e-05 B 33.373928 C 1 T 10.00 P 2 Prs -4.399772 Cnt -2317381 M 2.900796e-06 B 2.187160 C 1 T 10.00 P 3 Prs -3.850196 Cnt -2122031 M 2.900767e-06 B 2.187102 C 1 T 10.00 P 4 Prs -3.300011 Cnt -1926715 M 1.516154e-06 B -0.432091 C 1 T 10.00 P 5 Prs -2.750237 Cnt -1731705 M 3.562588e-05 B 51.719906 C 1 T 10.00 P 6 Prs -2.199967 Cnt -1537176 M 2.016204e-06 B 0.851326 C 1 T 10.00 P 7 Prs -1.650162 Cnt -1342878 M 2.529015e-06 B 1.487554 C 1 T 10.00 P 8 Prs -1.100179 Cnt -1148847 M 7.366149e-06 B 6.437026 C 1 T 10.00 P 9 Prs -0.549774 Cnt -955373 M 2.934924e-06 B 2.233797 C 1 T 10.00 P 10 Prs 0.000010 Cnt -762298 M 2.252551e-06 B 1.714509 C 1 T 10.00 P 11 Prs 0.549721 Cnt -569650 M 4.241780e-06 B 2.743043 C 1 T 10.00 P 12 Prs 1.099730 Cnt -377547 M 2.956774e-06 B 2.245238 C 1 T 10.00 P 13 Prs 1.649785 Cnt -185993 M 2.967257e-06 B 2.247349 C 1 T 10.00 P 14 Prs 2.199822 Cnt 5008 M 2.973132e-06 B 2.247442 C 1 T 10.00 P 15 Prs 2.749840 Cnt 195437 M 2.987920e-06 B 2.244944 C 1 T 10.00 P 16 Prs 3.299723 Cnt 336962 M 3.228561e-06 B 2.159962 C 1 T 10.00 P 17 Prs 3.849692 Cnt 574169 M 3.610811e-06 B 1.959847 C 1 T 10.00 P 18 Prs 4.399714 Cnt 762467 M 2.433504e-06 B 2.755404 C 1 T 10.00 P 19 Prs 4.949729 Cnt 950011 M 3.032422e-06 B 2.215327 C 1 T 10.00 P 20 Prs 5.499724 Cnt 1136719 M 3.032422e-06 B 2.215327 >

MAX-MIN CALIBRATION TABLE VALUES

DESCRIPTION	Returns the maximum and minimum calibration table values. This is used in determining overloads for scanning.
SYNTAX	MXMN
RETURNS	A prompt when ready for the next command.
EXAMPLE	MXMN 1 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 2 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 3 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 4 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 58 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 59 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 60 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 61 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 62 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 63 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 64 Max 5.5003 Min -5.5003 Overload Max 8.2504 Overload Min -8.2505 >

READ CURRENT TEMPERATURE

DESCRIPTION	Returns the current temperatures.
SYNTAX	TREAD [<channel>] Where channel is a temperature channel number from 1 to 8. When omitted all temperature channels are returned.
RETURNS	Temperature
EXAMPLE	TREAD <i>Temperature on sensor 1 is 38.625000</i> <i>Temperature on sensor 2 is 38.125000</i> <i>Temperature on sensor 3 is 38.625000</i> <i>Temperature on sensor 4 is 38.312500</i> <i>Temperature on sensor 5 is 38.250000</i> <i>Temperature on sensor 6 is 38.562500</i> <i>Temperature on sensor 7 is 38.375000</i> <i>Temperature on sensor 8 is 38.187500</i> > - OR- TREAD 1 <i>Temperature on sensor 1 is 37.062500</i> >

SCAN

DESCRIPTION	<p>Starts scanning and places the MPS4264 into SCAN mode.</p> <p>The MPS4264 uses the current PTP time as the scan start time under two conditions:</p> <ol style="list-style-type: none"> 1) When the PTP system is disabled via the SET PTPEN 0 command. 2) When the current time has not been received from the PTP server since the MPS4264 was powered up. Under this condition the MPS4264 time will start at Jan 1 1970, 0:00:00.000 hours, or epoch 0 time. This most likely will happen when there is no grandmaster on the network. <p>If the MPS was enabled and then disabled after receiving a time set from the PTP time server, the MPS will keep time from the last clock setting from the PTP server.</p> <p>When these two cases are not met, the MPS uses the start time set via the SET SST <time> and SET SSD <date> setting in the PTP group.</p> <p>When the current time is more recent than the scan start time, the MPS4264 will start scanning on the next frame interval that would be scheduled from the scan start time. See “Appendix D: 1588 Synchronized Scan Start” on page 85 for a discussion on starting scan using the PTP</p>
SYNTAX	SCAN
RETURNS	The Scan Data
EXAMPLE	SCAN NOTE: Output is based on formatting. >

MULTIPLE UNIT SCAN

DESCRIPTION	When the master receives the MSCAN command it issues a SCAN command to the multicast address set in the MCAST variable. Each MPS member of this multicast cluster will start scanning provided the SSD and SST conditions are met.
SYNTAX	MSCAN
RETURNS	The Scan Data
EXAMPLE	MSCAN NOTE: Output is based on formatting. >

SOFTWARE FRAME TRIGGER

DESCRIPTION	When in scan mode and external trigger set causes the MPS4264 to sample a frame of data.
SYNTAX	TRIG or <TAB>
RETURNS	A frame of data
EXAMPLE	TRIG >

VALIDATION

DESCRIPTION	Starts a validation and places the MPS into VAL mode
SYNTAX	<p>VAL [<F or FM>] Where: F- Initiates a single temperature field validation, with no oven connection. Connects to the calibrator. FM- Initiates a single temperature manual field validation with no oven or calibrator connection. When no parameter is present, the unit performs an automated validation based on the calibration group settings</p>
RETURNS	Validation Data
EXAMPLE	VAL >

MULTIPLE UNIT VALIDATION

DESCRIPTION	When the master receives the MVAL command it starts the multiple device validation sequence described in the section on multiple device calibration and validation.
SYNTAX	<p>MVAL [<F or FM>] Where: F- Initiates a single temperature multi-unit field validation, with no oven connection. Connects to the calibrator FM- Initiates a single temperature multi-unit manual field validation with no oven or calibrator connection. When no parameter is present, the unit performs an automatic multi-unit validation based on the calibration group settings.</p>
RETURNS	Validation Data
EXAMPLE	MVAL >

ZERO OFFSET CALIBRATION

DESCRIPTION	Starts a zero calibration and places the MPS4000 into CALZ mode.
SYNTAX	CALZ
RETURNS	A prompt when ready for the next command.
EXAMPLE	CALZ >

SAVE COMMANDS**SAVE**

DESCRIPTION	Saves all or a particular configuration to flash.																				
SYNTAX	<p>SAVE [<configuration>] When no configuration is present, all configurations except C and T are saved. When a configuration is present, only that group of variables is saved. The following table shows the configuration that are valid and their associated files on flash:</p> <table border="1"> <tr><td>IP</td><td>ip.cfg</td></tr> <tr><td>ID</td><td>id.cfg</td></tr> <tr><td>S</td><td>scan.cfg</td></tr> <tr><td>PTP</td><td>ptp.cfg</td></tr> <tr><td>T</td><td>Insert.txt</td></tr> <tr><td>C</td><td>Cal.cfg</td></tr> <tr><td>M</td><td>hw.cfg</td></tr> <tr><td>FTP</td><td>ftp.cfg</td></tr> <tr><td>FC</td><td>fc.cfg</td></tr> <tr><td>O</td><td>oven.cfg</td></tr> </table>	IP	ip.cfg	ID	id.cfg	S	scan.cfg	PTP	ptp.cfg	T	Insert.txt	C	Cal.cfg	M	hw.cfg	FTP	ftp.cfg	FC	fc.cfg	O	oven.cfg
IP	ip.cfg																				
ID	id.cfg																				
S	scan.cfg																				
PTP	ptp.cfg																				
T	Insert.txt																				
C	Cal.cfg																				
M	hw.cfg																				
FTP	ftp.cfg																				
FC	fc.cfg																				
O	oven.cfg																				
RETURNS	Prompt when saving is complete and ready for the next command.																				
EXAMPLE	To save all of the settings to all of the files (excluding C and T): SAVE >																				
EXAMPLE	To save only the “IP” settings to the file ip.cfg: SAVE IP – Save only the IP settings to the file ip.cfg >																				
NOTES	The “SAVE C” option saves all of the calibration coefficients along with the settings of the LIST C variable. As such, it can take up to 1 minute to complete.																				

SCAN VARIABLES**LIST SCAN SETTINGS (LIST S)**

DESCRIPTION	Returns all of the scan settings
SYNTAX	LIST S
RETURNS	All of the scan setting followed by the prompt.
EXAMPLE	LIST S <i>SET RATE 5.0000</i> <i>SET FPS 0</i> <i>SET UNITS PSI 1.000000</i> <i>SET FORMAT T F,F B,B B</i> <i>SET TRIG 0</i> <i>SET ENFTP 0</i> <i>SET OPTIONS 0 0 16</i> >

SET RATE (RATE)

DESCRIPTION	Set the scan rate in samples/channel/second (Hz)
SYNTAX	SET RATE <rate in HZ> Where rate is the number of samples per seconds for all 64 channels
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET RATE 5 >
DEFAULT	5
RANGE	.25 to 850

SET FRAMES PER SCAN (FPS)

DESCRIPTION	Set the number of frames in a scan
SYNTAX	SET FPS <number of frames> Where <number of frames> is the number of frames in a scan. When set to 0 the scan continues until the STOP command is issued.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET FPS 100 >
DEFAULT	0
RANGE	0 to 4294967295

SET SCAN UNITS (UNITS)

DESCRIPTION	Set the scan units
SYNTAX	<p>SET UNITS <units> [user units] Where units is one of the following:</p> <p>PSI ATM BAR CMHG CMH2O DECIBAR FTH2O GCM2 INHG INH2O KNM2 KGM2 KGCM2 KPA KIPIN2 MPA MBAR MH2O MMHG NM2 NCM2 OZIN2 OZFT2 PA PSF TORR USER RAW</p> <p>When USER units are selected, the user units must be supplied.</p>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET UNITS PSI >
EXAMPLE	SET UNITS USER 1.50 >
EXAMPLE	SET UNITS RAW >
DEFAULT	PSI

SET SCANNING DATA FORMAT (FORMAT)

DESCRIPTION	Set the format of the scanned data for each of the destinations: Telnet, FTP and Binary Server.
SYNTAX	<p>SET FORMAT <destination> <format code>[,<destination> <format code>,<destination> <format code>]</p> <p>Where:</p> <p><destination> is:</p> <ul style="list-style-type: none"> T-Telnet (only A, F and C format code apply) F-FTP (only A, C, B, and S format codes apply) B-Binary server (only B,L and S format codes apply) <p><format code> is:</p> <ul style="list-style-type: none"> A-ASCII F-Formatted ASCII C-CSV format B-Binary format (.dat) See binary format definition below. L-LabVIEW binary format. See binary format for LabVIEW definition below. S-Statistical data (.dat) See statistical binary format definition below.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET FORMAT T F, F B, B B >
DEFAULT	T F, F B, B B

SET TRIGGER (TRIG)

DESCRIPTION	Sets the frame trigger
SYNTAX	<p>SET TRIG <trig></p> <p>When trig is set to 0 frame trigger is internal.</p> <p>When trig is set to 1 frame trigger is external or triggered with the TRIG command or <TAB> key.</p> <p>When trig is set to 2, external scan trigger is configured.</p> <p>When trig is set to 3 auto scan is enabled (see "Auto Scan Option" on page 24)</p>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET TRIG 0 >
DEFAULT	0
RANGE	0, 1, 2 or 3
NOTES	Setting TRIG to 2 requires a SAVE S command followed by a reboot in order for the change to take effect.

SET ENABLE FTP (ENFTP)

DESCRIPTION	Enables or disables sending data to the users FTP server.
SYNTAX	SET ENFTP <0-Disable, 1-Enable>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET ENFTP 1 >
DEFAULT	0
RANGE	0 or 1

SET OPTIONS (OPTIONS)

DESCRIPTION	Set the subset size of the real time data analysis sample
SYNTAX	<p>SET OPTIONS <fast scan> <read mode> <statistical scan subset size></p> <p>Where <fast scan>: 0 – Normal scan 1 through 4 – Fast scan starting channel (see section on “Fast Scan Option (850Hz - 2500Hz)” on page 24)</p> <p><read mode>: 0 – Normal statistical scan (see section on “Statistical Scan Option” on page 24) 1 – Read mode in statistical scan</p> <p><statistical scan subset size>: Size is set to any value between 2 and 256. This value is the number of samples the statistical engine uses to calculate the statistical parameters.</p>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET OPTIONS 0 0 16 >
DEFAULT	0 0 16
RANGE	0 through 4, 0 or 1, 2 to 256

UDP VARIABLES

The UDP settings control the target and enable of scan data out the UDP port.

LIST UDP SETTINGS (LIST UDP)

DESCRIPTION	Returns all of the UDP settings
SYNTAX	LIST UDP
RETURNS	All of the setting followed by the prompt.
EXAMPLE	LIST UDP <i>SET ENUDP 0</i> <i>SET IPUDP 224.0.1.2</i> >

SET UDP ENABLE (ENUDP)

DESCRIPTION	Sets the UDP Enable
SYNTAX	SET ENUDP <0 or 1>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET ENUDP 1 >
DEFAULT	0
RANGE	0 or 1

SET THE UDP TARGET ADDRESS AND PORT (IPUDP)

DESCRIPTION	Sets the UDP target address and the port. After setting this value a SAVE must be issued and the MPS4000 must be power cycled or rebooted for the changes to take effect.
SYNTAX	SET IPUDP <udpaddr> <port> Where: <udpaddr> - Is the target address for the UDP scan data. <port> - Is the target port.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET IPUDP 224.0.1.2 23 >
DEFAULT	0.0.0.0

IDENTIFICATION VARIABLES

The identification settings control the serial number and the nominal pressure range of the MPS4000.

LIST IDENTIFICATION SETTINGS (LIST ID)

DESCRIPTION	Returns all of the identification settings
SYNTAX	LIST ID
RETURNS	All of the setting followed by the prompt.
EXAMPLE	LIST ID <i>SET SN 222</i> <i>SET NPR 5.0000 -5.0000 5.0000 -5.0000</i> <i>SET MCAST 224.0.1.2</i> <i>></i>

SET SERIAL NUMBER (SN)

DESCRIPTION	Sets the serial number of the MPS
SYNTAX	SET SN <sn> Where: <sn> is the serial number ranging from 0 to 32767.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET SN 123 <i>></i>
DEFAULT	100
RANGE	0 to 32767

SET NOMINAL PRESSURE RANGE (NPR)

DESCRIPTION	Sets the nominal pressure ranges for the MPS
SYNTAX	SET NPR <npr1p> <npr1n> <npr2p> <npr2n> Where: <npr1p> is the maximum nominal pressure for channels 1 to 32. <npr1n> is the minimum nominal pressure for channels 1 to 32. <npr2p> is the maximum nominal pressure for channels 33 to 64. <npr2n> is the minimum nominal pressure for channels 33 to 64.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET NPR 15.0 -15.0 15,0 -15,0 <i>></i>
DEFAULT	15.0 15.0 15.0 15.0

SET THE MULTIPLE DEVICE MULTICAST ADDRESS (MCAST)

DESCRIPTION	Sets the multicast IP address used for multiple unit calibration and validation. This address must be the same for all units to be calibrated as a group.
SYNTAX	SET MCAST <ipadd> Where: <ipadd> - Is the multicast address in the range of 224.0.0.0 to 239.255.255.255
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET MCAST 224.0.1.2 <i>></i>
DEFAULT	224.1.1.11
NOTES	A save and power cycle is required for changes to take effect.

IP VARIABLES

The IP settings control the IP address and associated network identification settings of the MPS4000. The IP configuration is accessed by the command LIST IP.

LIST IP SETTINGS (LIST IP)

DESCRIPTION	Returns all of the IP settings
SYNTAX	LIST IP
RETURNS	All of the IP setting followed by the prompt.
EXAMPLE	LIST IP <i>SET IPADD 10.0.1.222</i> <i>SET SUBNET 255.255.0.0</i> <i>SET MAC 0.96.93.90.0.100</i> <i>SET LOGIN Scanivalve</i> <i>SET PASSWORD Scanner</i> <i>SET LOGIN1 Scanivalve1</i> <i>SET PASSWORD1 Scanner1</i> <i>SET ALLOWANON 1</i> <i>SET APP Mps4000.hex</i> <i>SET GW 0.0.0.0</i> <i>SET NAME MPS0000</i> >

SET IP ADDRESS OF THE MPS (IPADD)

DESCRIPTION	Sets the IP address of the MPS4000
SYNTAX	SET IP <ipadd> Where <ipadd> is any valid IP address in the form nnn.nnn.nnn.nnn When <ipadd> is set to 0.0.0.0 then DHCP is used. See the section on DHCP.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET IPADD 10.0.1.222 >
DEFAULT	191.30.105.100
RANGE	Any valid IP address

SET SUBNET OF THE MPS (SUBNET)

DESCRIPTION	Sets the subnet mask for the MPS4000
SYNTAX	SET SUBNET <mask> Where <mask> is in the form nnn.nnn.nnn.nnn
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET SUBNET 255.255.0.0 >
DEFAULT	255.255.255.0
RANGE	Any valid subnet mask

SET LOGIN NAME (LOGIN)

DESCRIPTION	Sets the accepted login user name for the FTP server in the MPS4000
SYNTAX	SET LOGIN <user name> Where <user name> is any ASCII string with out spaces up to 64 characters long.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET LOGIN SCANIVALVE >
DEFAULT	Scanivalve

SET PASSWORD FOR LOGIN NAME (PASSWORD)

DESCRIPTION	Sets the password for the accepted user names in LOGIN for the FTP server in the MPS 4000
SYNTAX	SET PASSWORD <password> Where <password> is any ASCII string with out spaces up to 64 characters long.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET PASSWORD scanner >
DEFAULT	Scanner

SET LOGIN NAME 1 (LOGIN1)

DESCRIPTION	Sets an additional accepted login user name for the FTP server in the MPS4000
SYNTAX	SET LOGIN1 <user name> Where <user name> is any ASCII string with out spaces up to 64 characters long.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET LOGIN1 SCANIVALVE1 >
DEFAULT	Scanivalve1

SET PASSWORD FOR LOGIN NAME 1 (PASSWORD1)

DESCRIPTION	Sets the password for the accepted user names in LOGIN1 for the FTP server in the MPS 4000
SYNTAX	SET PASSWORD1 <password> Where <password> is any ASCII string with out spaces up to 64 characters long.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET PASSWORD1 scanner1 >
DEFAULT	Scanner1

SET ALLOW ANONYMOUS FTP LOGIN (ALLOWANON)

DESCRIPTION	Allows or prevents the FTP server to allow anonymous logins. When set to 1 the server accepts anonymous login when set to 0 the user must login with one of the two login name and password combination in LOGIN and LOGIN1.
SYNTAX	SET ALLOWANON <1 or 0>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET ALLOWANON 1 >
DEFAULT	1
0 or 1	0 or 1

SET APPLICATION FILE NAME (APP)

DESCRIPTION	The file name of the application to run. This is the file name that is used when automatically running the application from the boot loader. It is also the file name used when using the RUN command. If this file is not found, an error is returned. File name is case sensitive.
SYNTAX	SET APP <file name> Where <file name> is the file name to run on start up. It is case sensitive.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET APP Mps4000.hex >
DEFAULT	Mps4000.hex

SET GATEWAY (GW)

DESCRIPTION	This is the IP address of the network gateway used for network communication. If not present 0.0.0.0 must be used.
SYNTAX	SET GW <ip address> Where <ipaddr> is any valid IP address in the form nnn.nnn.nnn.nnn
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET GW 10.0.2.123 >
DEFAULT	0.0.0.0

SET DNS NAME (NAME)

DESCRIPTION	This sets the name sent to the network DNS server to identify this MPS4000. When DNS name resolution is used.
SYNTAX	SET NAME <DNS name> Where <DNS name> is any ASCII string up to 64 characters. Spaces are not allowed.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET NAME MPS0222 >
DEFAULT	MPS0000

PRECISION TIME PROTOCOL VARIABLES (LIST PTP)

The precision time protocol settings control the PTP 1588 settings of the MPS4000. The PTP configuration is accessed by the command LIST PTP.

LIST PRECISION TIME PROTOCOL SETTINGS (LIST PTP)

DESCRIPTION	Returns all of the PTP settings
SYNTAX	LIST PTP
RETURNS	All of the PTP setting followed by the prompt.
EXAMPLE	<pre>LIST PTP SET PTPEN 1 SET STAT 0 SET SST 0:0:0.000000 SET SSD 1971/1/1 SET UTCOFFSET -7:00:00 ></pre>

SET PTP ENABLE (PTPEN)

DESCRIPTION	Enables the precision time protocol engine in the MPS4000. When PTPEN is set to 2, the MPS will serve as a PTP master
SYNTAX	SET PTPEN <0-disable, 1-enable, 2-enable as master>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	<pre>SET PTPEN 1 ></pre>
DEFAULT	0
RANGE	0, 1 or 2

SET PTP STATISTICAL OUTPUT (STAT)

DESCRIPTION	Sets the destination of PTP statistical output
SYNTAX	<pre>SET STAT <level code> Where <level code> is: 0 – off 1 – serial 2 – network</pre> <p>The statistical data is output in the following format: <ofm>, <msd>, <smd>, <mpd> ofm - Offset From Master. Time that slave and master differ in units of microseconds msd - Master to Slave Delay. Filtered packet delivery time from master to slave smd - Slave to Master Delay. Filtered packet delivery time from slave to master mpd - Mean Path Delay. Filtered mean path delay time</p>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	<pre>SET STAT 0 ></pre>
DEFAULT	0
RANGE	0 to 2

SET SCAN START TIME (SST)

DESCRIPTION	Sets the time to start scanning as referenced to the grandmaster. The MPS4000 will use the current PTP time as the scan start time under two conditions: 1) When the PTP system is disabled via the SET PTPEN 0 command. 2) When the current time has not been received from the PTP server since the MPS4000 was powered up. Under this condition the MPS4000 time will start at Jan 1 1970, 0:00:00.000 hours, or epoch 0 time. This most likely will happen when there is no PTP server on the network. If the MPS4000 was enabled and then disabled after receiving a time set from the PTP time server, the MPS4000 will keep time from the last clock setting from the PTP server. When these two cases are not met, the MPS4000 uses the start time set via the SET SST <time> and SET SSD <date> setting in the PTP group.
SYNTAX	SET SST <hh:mm:ss.fffff> Where: hh-Hours, mm- Minutes, ss-Seconds and ffff is fractions of a second.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET SST 13:00:00.000 >
DEFAULT	0:0:0.0000
RANGE	Any valid time

SET SCAN START DATE (SSD)

DESCRIPTION	Sets the date to start scanning as referenced to the grandmaster.
SYNTAX	SET SSD <yyyy/mm/dd > Where: yyyy-Year, mm- Month, dd- Day
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET SSD 2016/08/10 >
DEFAULT	1971/1/1
RANGE	Any valid date

SET UTC OFFSET (UTCOFFSET)

DESCRIPTION	Sets the offset to be added or subtracted from network time.
SYNTAX	SET UTCOFFSET hh:mm:ss Where: hh-Hours from -12 to 12 mm- Minutes from 0 to 59 ss – Seconds from 0 to 59
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET UTCOFFSET -9:0:0 >
DEFAULT	0:0:0
RANGE	Any valid time

CONVERSION TABLE VARIABLES

The conversion table is responsible for converting raw pressure counts and temperature into accurate EU pressure.

LIST CONVERSION TABLE SETTINGS (LIST T)

DESCRIPTION	Returns the selected conversion table settings
SYNTAX	LIST T [<start temp>] [<end temp>] [<channel>]
RETURNS	All of the selected conversion table settings in the format: SET PT <temp> <channel> <pressure> <counts> <cal time> Where: <temp> - is the temperature of the calibration point <channel> - is the channel form 1 to 64 <pressure> - is the pressure of the calibration point <counts> - is the A/D counts at the calibration pressure <cal time> - is the PTP time of the calibration in seconds from 1900.
EXAMPLE	LIST T 35 40 64 SET PT 38.25 64 -5.499972 2609756 370714637 SET PT 38.26 64 -4.124957 1714679 370714637 SET PT 38.26 64 -2.750003 801797 370714637 SET PT 38.27 64 -1.375020 -127217 370714637 SET PT 38.29 64 -0.000031 -1070892 370714637 SET PT 38.27 64 1.374991 -2026530 370714637 SET PT 38.27 64 2.750021 -2990477 370714637 SET PT 38.27 64 4.125011 -3961499 370714637 SET PT 38.28 64 5.500008 -4936479 370714637 >

SET CONVERSION TABLE POINT (PT)

DESCRIPTION	Sets a conversion table point
SYNTAX	SET PT <temp> <channel> <pressure> <counts> [<cal time>] Where: <temp> - is the temperature of the calibration point <channel> - is the channel form 1 to 64 <pressure> - is the pressure of the calibration point <counts> - is the A/D counts at the calibration pressure <cal time> - is the PTP time of the calibration in seconds from 1900.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET PT 38.28 64 5.500008 -4936479 >

CALIBRATION VARIABLES (LIST C)

The calibration settings control the calibration and validation of the MPS4000. Reference Scanivalve Document: MPS Field Val if performing a single temperature validation with an SPC4000.

LIST CALIBRATION SETTINGS (LIST C)

DESCRIPTION	Returns all of the calibration settings
SYNTAX	LIST C
RETURNS	All of the calibration settings followed by the prompt.
EXAMPLE	LIST C <i>SET NUMPTS 0 9 0 9</i> <i>SET MIN 40.000000 -5.500000 5.000000 -5.000000</i> <i>SET MAX 60.000000 5.500000 55.000000 5.000000</i> <i>SET DELAY 5 5</i> <i>SET IPOVEN 0.0.0.0 0</i> <i>SET IPCAL 0.0.0.0 0 1</i> <i>SET CALAVG 1 1</i> <i>SET VALZO 0</i> <i>SET FCAL 3 -5.00 5.00</i> <i>SET FVAL 3 -5.00 5.00</i> <i>></i>

SET NUMBER OF POINTS IN CONVERSION TABLE (NUMPTS)

DESCRIPTION	Sets the number calibration and validation temperature and pressure points. This setting works with the MIN and MAX setting to the calibration and validation points.
SYNTAX	SET NUMPTS <tcal> <pcal> <tval> <pval> Where: <tcal> - Number of temperature points for a calibration <pcal> - Number of pressure points for a calibration <tval> - Number of temperature points for a validation <pval> - Number of pressure points for a validation. These values are not range checked.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET NUMPTS 5 9 5 9 >
DEFAULT	5 9 5 9
RANGE	<0 to 15> <0 to 25> <0 to 15> <0 to 25>

SET MINIMUM CALIBRATION AND VALIDATION VALUES (MIN)

DESCRIPTION	Sets the minimum calibration and validation temperature and pressure.
SYNTAX	SET MIN <tcal> <pval> <tval> <pval> Where: <tcal> - Minimum temperature value for a calibration. <pval> - Minimum pressure value for a calibration <tval> - Minimum temperature value for a validation <pval> - Minimum pressure value for a validation.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET MIN -5.5 -5.5 -5.0 -5.0 >
DEFAULT	0 -15 0 -15
RANGE	These values are not range checked.

SET MAXIMUM CALIBRATION AND VALIDATION VALUES (MAX)

DESCRIPTION	Sets the maximum calibration and validation temperature and pressure.
SYNTAX	SET MAX <tcal> <pval> <tval> <pval> Where: <tcal> - Maximum temperature value for a calibration. <pval> - Maximum pressure value for a calibration <tval> - Maximum temperature value for a validation <pval> - Maximum pressure value for a validation.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET MAX 5.5 5.5 5.0 5.0 >
DEFAULT	70 15 70 15
RANGE	These values are not range checked.

SET OVEN AND CALIBRATOR DELAYS (DELAY)

DESCRIPTION	Sets the delay in minutes and seconds for the oven and calibrator, respectively.
SYNTAX	SET DELAY <iodelay> <odelay> <cdelay> Where: <iodelay> - Initial oven delay in minutes <odelay> - Oven delay in minutes. Setting this value to zero enables a stable query to the MPS (MPS must be +/-0.125°C for 10 minutes, 120 minute timeout). <cdelay> - Calibrator delay in seconds. Setting this value to zero uses the calibrator's STABLE command.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET DELAY 60 15 >
DEFAULT	60 45
RANGE	These values are not range checked.

SET THE CALIBRATOR IP ADDRESS, PORT AND CALIBRATOR NUMBER (IPCAL)

DESCRIPTION	Sets the IP address and port for the calibrator. When 0.0.0.0 0 0 is set a simulated calibrator for debug is used. This is only applicable when using a Scanivalve SPC4000 calibrator.
SYNTAX	SET IPCAL <ipadd> <port> <calibrator number> Where: <ipadd> - Calibrator IP address <port> - Calibrator port <calibrator number> - Calibrator number
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET IPCAL 10.0.0.122 23 1 >
DEFAULT	10.0.0.61 23 1
RANGE	These values are not range checked.

SET THE CALIBRATION AND VALIDATION AVERAGE (CALAVG)

DESCRIPTION	Sets the number of samples taken from the MPS for a calibration and validation.
SYNTAX	SET CALAVG <calavg> <valavg> Where: <calavg> - Is the number of sample taken for a calibration. Valid values range between 1 and 32000. <valavg> - Is the number of sample taken for a validation. Valid values range between 1 and 32000.
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET CALAVG 16 16 >
DEFAULT	1
RANGE	<1 to 32000> <1 to 32000>

SET ENABLE CALZ BEFORE VALIDATION (VALZO)

DESCRIPTION	When set a CALZ is performed at the start of each temperature during a validation.
SYNTAX	SET VALZO <0 or 1>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET VALZO 1 >
DEFAULT	0
RANGE	0 or 1

SET FIELD CALIBRATION POINTS (FCAL)

DESCRIPTION	Sets the number of points and the min/max pressure for field calibration. Points will be evenly spaced between the minimum and maximum and will include a 0.00 point.
SYNTAX	SET FCAL <number of points> <minimum pressure> <maximum pressure>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET FCAL 5 -5.00 5.00 >
DEFAULT	0
RANGE	These values are not range checked

SET FIELD VALIDATION POINTS (FVAL)

DESCRIPTION	Sets the number of points and the min/max pressure for field validation. Points will be evenly spaced between the minimum and maximum and will include a 0.00 point.
SYNTAX	SET FVAL <number of points> <minimum pressure> <maximum pressure>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET FVAL 3 -5.00 5.00 >
DEFAULT	0
RANGE	These values are not range checked

OVEN VARIABLES (LIST O)

The oven settings control the thermal chamber when performing a multi-temperature calibration or validation of the MPS4000. Reference Scanivalve Document: MPS Field Val if performing a single temperature validation with an SPC4000.

LIST OVEN SETTINGS (LIST O)

DESCRIPTION	Returns all of the oven settings
SYNTAX	LIST C
RETURNS	All of the oven settings followed by the prompt.
EXAMPLE	LIST O <i>SET IPOVEN 0.0.0.0 0</i> <i>SET STARTOVEN 0</i> <i>SET STOPOVEN 0</i> <i>SET TEMPOVEN 0</i> <i>></i>

SET THE OVEN IP ADDRESS, PORT AND CALIBRATOR NUMBER (IPOVEN)

DESCRIPTION	Sets the IP address and port for the oven When 0.0.0.0 0 0 is set no oven is used. This is only applicable when using performing multi-temperature calibrations or validations.
SYNTAX	SET IPOVEN <ipadd> <port> Where: <ipadd> - Calibrator IP address <port> - Calibrator port
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET IPOVEN 10.0.0.122 23 1 <i>></i>
DEFAULT	0.0.0.0 0
RANGE	These values are not range checked.

SET THE OVEN START COMMAND (STARTOVEN)

DESCRIPTION	Sets the command to send the oven to start
SYNTAX	SET STARTOVEN [<start oven command>] NOTE: If no command is entered, a CR is sent to the oven
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET STARTOVEN RUNM <i>></i>
DEFAULT	0
RANGE	These values are not range checked.

SET THE OVEN STOP COMMAND (STOPOVEN)

DESCRIPTION	Sets the command to stop the oven.
SYNTAX	SET STOPOVEN [<stop oven command> NOTE: If no command is entered, a CR is sent to the oven]
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET STOPOVEN STOP >
DEFAULT	0
RANGE	These values are not range checked.

SET THE TEMPERATURE COMMAND (TEMPOVEN)

DESCRIPTION	Sets the command to send the oven to a temperature
SYNTAX	SET TEMPOVEN [<set temperature command> NOTE: If no command is entered, a CR is sent to the oven]
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET TEMPOVEN SETP1, >
DEFAULT	0
RANGE	These values are not range checked.

MISCELLANEOUS VARIABLES (LIST M)

This is a collection of settings that are used for debugging or other setting.

LIST MISCELLANEOUS SETTINGS (LIST M)

DESCRIPTION	Returns all of the miscellaneous settings
SYNTAX	LIST M
RETURNS	All of the miscellaneous setting followed by the prompt.
EXAMPLE	LIST M SET SIM 0 SET ECHO 0 SET XITE 2 0 1 SET SVRSEL 1 >

SET SIMULATION (SIM)

DESCRIPTION	Sets debugging and simulation modes.
SYNTAX	SET SIM <code in hex> Where <code in hex>: Bit 2 when set shows PTP calculated start differentials Bit 0 when set uses a prototype board
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET SIM 4 >
DEFAULT	0
RANGE	0 to 65535

SET ECHO (ECHO)

DESCRIPTION	Sets echo for serial and Telnet communication.
SYNTAX	SET ECHO <0- no echo, 1-echo >
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET ECHO 1 >
DEFAULT	0
RANGE	0 or 1

SET EXCITATION (XITE)

DESCRIPTION	Sets excitation mode
SYNTAX	SET XITE <0-neg only, 1-pos only, 2-reg, 3-reversed polarity> <man A/D rate 0 for auto> <A/D Idle enabled>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET XITE 2 0 1 >
DEFAULT	2 0 1
RANGE	<0,1 or 2> <0 to 9 or F> <0, 1>

SET SERVER SELECT (SVRSEL)

DESCRIPTION	Sets the type of network server. See section on socket allocation. After changing any of these settings a SAVE command must be issued and the device must be restarted.
SYNTAX	SET SVRSEL <1, 2, or 3> Where: 1 – Default Multicast server 2 – Binary server 3 – UDP scan data out
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET SVRSEL 1 >
DEFAULT	1
RANGE	1, 2 or 3

FTP VARIABLES (LIST FTP)

The FTP settings are used to control the MPS4000 FTP client connection to the FTP server.

LIST FTP SETTINGS (LIST FTP)

DESCRIPTION	Returns all of the FTP settings
SYNTAX	LIST FTP
RETURNS	All of the FTP setting followed by the prompt.
EXAMPLE	LIST FTP <i>SET USERFTP Scanivalve</i> <i>SET PASSFTP password</i> <i>SET PATHFTP C:/FTPDIR</i> <i>SET IPFTP 10.0.1.222</i> <i>SET FILEFTP Scan.txt</i> >

SET USER NAME (USERFTP)

DESCRIPTION	Sets the name of the user login on the FTP server or NAS
SYNTAX	SET USERFTP <user name>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET USERFTP Scanivalve >
DEFAULT	Admin

SET FTP PASSWORD (PASSFTP)

DESCRIPTION	Sets the name of the user password for login on the FTP server or NAS
SYNTAX	SET PASSFTP <password>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET PASSFTP password >
DEFAULT	Password

SET FTP PATH (PATHFTP)

DESCRIPTION	Set the directory path of the scan data file on the user's computer.
SYNTAX	SET PATHFTP <directory path>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET PATHFTP /disk1/share >
DEFAULT	/disk1/share

SET FTP IP ADDRESS (IPFTP)

DESCRIPTION	Set the IP address of the users FTP server.
SYNTAX	SET IPFTP <IP Address>
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET PATHFTP 10.0.1.222 >
DEFAULT	10.0.0.1

SET FTP FILE NAME (FILEFTP)

DESCRIPTION	Set the base name for the scan data on the user's computer. The actual name will have the time and date appended to the file. It will have a file extension based on the data format.
SYNTAX	SET FILEFTP scan
RETURNS	Prompt when ready to accept a new command
EXAMPLE	SET FILEFTP Scan_117_ >
DEFAULT	SCAN

[INTENTIONALLY LEFT BLANK]

SECTION 6: MAINTENANCE

FIRMWARE INSTALLATION & CONFIGURATION

FILE UPLOAD (WINDOWS XP, 7, 10)

NOTE: Some firewalls may block file transfers. We recommend that firewalls be shut down for the procedures in this section.

1. Connect a host computer to the Ethernet port.
2. Locate the new or updated firmware (MpsXXX.TAR)
3. Open My Computer
4. In the Address Bar, type:
ftp://<IPAddress><Enter>
Where: <IPAddress> is the IP Address of the MPS4264.
5. The host computer should connect, If the connection is successful, the contents of the MPS4264 flash storage will be displayed in a folder format. Find and delete the old MPSXXX.tar file.

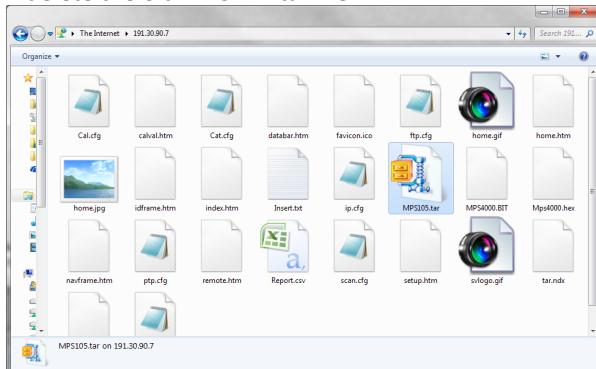


FIGURE 6-1: FLASH FILE CONTENTS

6. Open Windows Explorer in another window. Find the directory where the new MpsXXX.TAR file is stored.
- Highlight the file and

Click: Edit


Click: Copy

7. Select the MPS4264 window and,

Click: Edit

Click: Paste

The file will be copied to the MPS4264. A progress window will provide progress information.

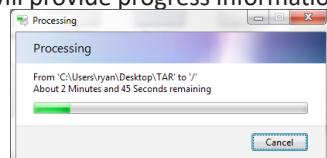

8. Once the file has been copied, connect to the MPS.
- To "unpack" the .TAR file, send the command:

```
TAR <FILENAME>
```

The MPS will display the status as it unpacks the file.

```
ScanTel
File Disconnect Configuration Font Size Display Help
>tar MPS4264.tar
Starting to de-tar files from fileMPS4264.tar
Extracting file MPS4000.BIT 340701 2016/05/18 18:29:30
Extracting file MPS4000.hex 695026 2016/06/13 20:15:41:24
Extracting file calval.htm 19639 2015/04/17 15:05:33
Extracting file database.htm 18118 2016/01/12 18:56:09
Extracting file favicon.ico 1150 2014/03/14 15:19:55
Extracting file home.gif 81414 1980/01/01 19:00:00
Extracting file home.jpg 16814 2014/05/12 19:14:03
Extracting file home.htm 761 2015/05/22 22:51:24
Extracting file idframe.htm 3933 2015/06/17 17:08:13
Extracting file index.htm 2087 2015/05/27 14:47:55
Extracting file navframe.htm 2994 2015/05/27 16:36:18
Extracting file remote.htm 8141 2015/04/18 19:48:57
Extracting file setup.htm 10312 2016/03/08 16:40:52
Extracting file svLogo.gif 2472 2015/06/01 17:03:21
Extracting file term.htm 3101 2015/06/03 19:33:35
Done extracting files from MPS4264.tar
Creating TAR.NDX file
Done creating TAR.NDX file
IP Address 191.30.90.7 | Connect Status Connected | File Capture Off
```

FIGURE 6-2: UNPACKING TAR FILE

Once it is completed, verify the files were unpacked correctly by sending the command:

```
TARCHK
```

The MPS will display a confirmation that the files were unpacked successfully.

```
ScanTel
File Disconnect Configuration Font Size Display Help
>tarchk
Starting TARCHK
Testing file MPS4000.BIT
Testing file MPS4000.hex
Testing file calval.htm
Testing file FavIcon.ico
Testing file home.gif
Testing file home.jpg
Testing file home.htm
Testing file idframe.htm
Testing file navframe.htm
Testing file remote.htm
Testing file setup.htm
Testing file svLogo.gif
Testing file term.htm
PMS4264.tar: 0 errors
IP Address 191.30.90.7 | Connect Status Connected | File Capture Off
```

FIGURE 6-2: UNPACKING TAR FILE

(See "TAR File Unpack" on page 43 and "Firmware Files Check" on page 41 for more information on the TAR and TARCHK commands.)

9. Cycle the AC power or execute a Reboot command to complete the process.

CALIBRATION COEFFICIENT INSTALLATION (WINDOWS XP, 7, 10)

NOTE: Some firewalls may block file transfers. We recommend that firewalls be shut down for the procedures in this section.

1. Connect a host computer to the Ethernet port.
 2. Open My Computer
 3. In the Address Bar, type:
ftp://<IPAddress><Enter>
- Where: <IPAddress> is the IP Address of the MPS4264.
4. The host computer should connect. If the connection is successful, the contents of the MPS4264 flash storage will be displayed in a folder format.

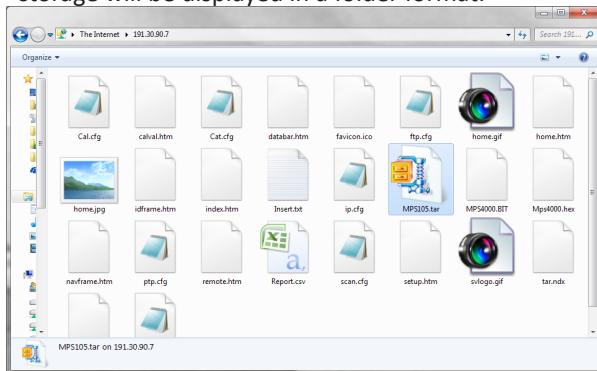


FIGURE 6-4: FLASH FILE CONTENTS

5. Open Windows Explorer in another window. Find the directory where the new calibration coefficients are stored. MPS4264 calibration coefficients are a text file called "CAL.CFG". Highlight the files to be copied to the MPS4264 and

Click: Edit


Click: Copy

6. Select the MPS4264 window and,

Click: Edit

Click: Paste

The files will be copied to the MPS4264, replacing any existing CAL.CFG file that may have existed. A progress window will provide progress information.

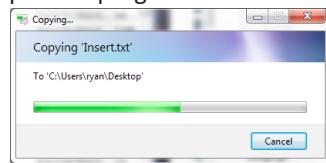

7. When the copy function is complete, the MPS4264 window will show the new files. The time and date of the new file can be used to confirm that the new files were uploaded successfully.

8. Cycle the AC power or execute a Reboot command to complete the process.

DIP SWITCHES & BOOTLOADER MODE

DIP SWITCHES

The MPS4264 has two DIP switches on the Power/Ethernet board. The configuration of the first DIP switch enables or disables detailed debug output from the unit. The configuration of the second DIP switch determines the module's function regarding when the module boots into the "application" file and when it is held in "bootloader" mode. The two switches are:

SW1 - When this switch is on, the boot loader will run in the debug mode. Debug output is directed to the serial port. Default is off.

SW2 - When this switch is on, the application automatically boots on power up. When this switch is off, the MPS4264 remains in the bootloader mode. Default is on.

*NOTE: The "ON" and "OFF" positions are the opposite of the "ON" label on the switches. Figure 6-6 shows the DIP switches with the "correct" position labels.

WHEN TO OPERATE IN BOOTLOADER MODE

If the MPS will not boot, or appears to have lost communications with a Host, the system may be operated in the Bootloader Mode. This is a "safe mode" of sorts where the MPS stays in the "boot loader." In this mode of operation, the application will not start. All available commands may be entered by the Ethernet or Configuration ports. Typically, the only common reason to use this function is if the application file has been corrupted and the MPS is failing to boot from it correctly.


HOW TO OPERATE IN BOOTLOADER MODE

1. Power the MPS down and disconnect the power and Ethernet cable. Ensure proper ESD considerations are taken.




CAUTION! ESD PROTECTION REQUIRED. The proper use of grounded work surfaces and personal wrist straps are required when coming into contact with exposed circuits to prevent static discharge from damaging sensitive electronic components.

2. Remove the two screws that mount the end cover in place. Carefully remove the end cover by sliding it vertically up. Care must be taken to not bend or move the power or Ethernet connector.

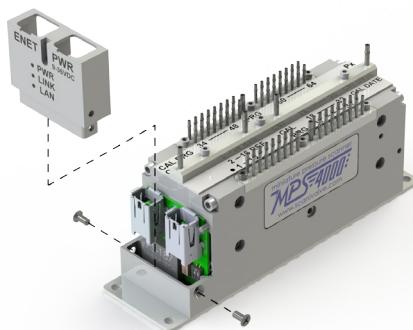**FIGURE 6-5: END COVER REMOVAL**

- With the end cover removed, the DIP switches are accessible. To hold the unit in the bootloader and prevent it from booting into the application file, switch the #2 DIP switch downwards to the “OFF” position. NOTE: the switch is labeled with an “ON” label but for this application it is backwards. “ON” is up and “OFF” is down.

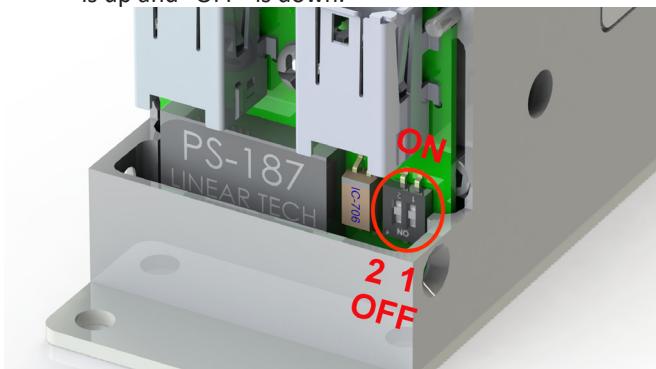**FIGURE 6-6: DIP SWITCHES**

- Reconnect the power and Ethernet cables. Re-apply power to the MPS4264. The application will not run, but the bootloader commands can be used to troubleshoot system problems.
- When troubleshooting is complete, remove power and reset switch #2 to the “ON” position. Reinstall the end cover and re-apply power.

FORMATTING THE FLASH DISK

Should the flash disk of the MPS need to be reformatted, it can be done fairly simply if the following procedure is carefully followed. For most, but not all times that an MPS will be reformatted a serial connection needs to be available to communicate with the unit. If you are able to currently establish an Ethernet connection with the unit, the reformatting is not the best option. In that case, simply FTP to the unit, manually delete the contents of the directory and reload them using the procedure outlined in “Firmware Installation & Configuration File Upload (Windows XP, 7, 10)” on page 76. However, if an Ethernet connection cannot be established with the unit it is probably the case

that the MPS application file or memory system has been corrupted and needs to be reformatted. The procedure for reformatting the MPS is as follows:

- Make sure the MPS is powered down and proper ESD considerations are taken.
- Following the procedure outlined in “How to operate in bootloader mode” on page 77, configure the DIP switch for bootloader mode with positioning switch #2 to “OFF.”
- Power the MPS on. If the IP address is known, you can attempt to connect to the unit via an Ethernet connection. If this is not successful, or if the IP address is not known you will need to connect using the serial port. Details on connecting to the serial port can be found in “Serial Communications” on page 14.
- Once a connection is established, send the command:


FDISK

The MPS will respond with “Type FDISKCONFIRM to confirm FDISK or STOP to escape” Follow the instructions and send the command:

FDISKCONFIRM

The MPS will respond indicating completion:

```
Start format
Erasing device 100
Done device 100
```

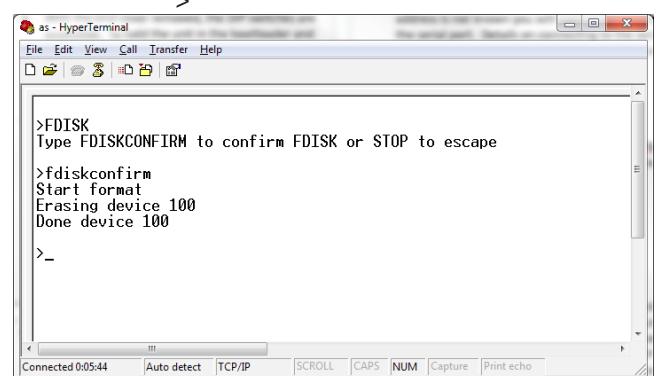**FIGURE 6-7: FDISK COMMAND**

5. At this point the flash disk has been reformatted and all contents have been deleted. All configuration settings, including the MPS's IP address have been deleted. Before the power is cycled the IP address needs to be set so Ethernet communications can be established. Send the command:


SET IPADD XXX.XXX.XXX.XXX

(where xxx.xxx.xxx.xxx is the desired IP address)

Followed by the command:

SAVEIP

SAVEIPCONFIRM

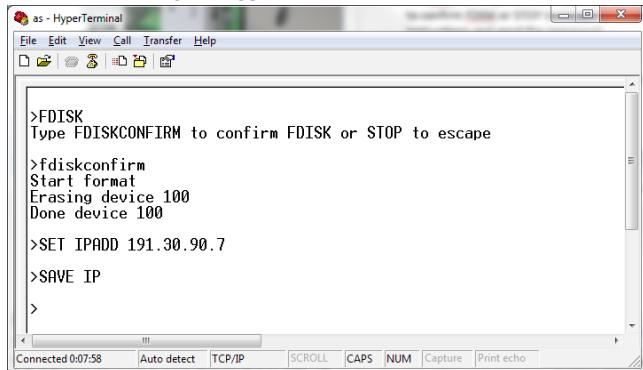

```

>FDISK
Type FDISKCONFIRM to confirm FDISK or STOP to escape
>fdiskconfirm
Start format
Erasing device 100
Done device 100
>SET IPADD 191.30.90.7
>SAVE IP
>

```

FIGURE 6-8: RESET IP ADDRESS

6. FTP to the MPS and drag-and-drop the two files:

Mps4000.hex

Mpsxxxx.tar

- using the procedure outlined in “Firmware Installation & Configuration File Upload (Windows XP, 7, 10)” on page 76.
7. Power the MPS down, and reconfigure the DIP switches for normal operation with switch #1 “OFF” and switch #2 “ON.”
 8. Power the MPS back up and connect via the Ethernet connection as normal using the IP Address set in step # 5.
 9. Send the command:

TAR MPSxxx.TAR

(using the correct file name of the .TAR file)

- This “unpacks” the .TAR file and sets up the file structure of the MPS.
10. At this point, the application files and web server files have been restored. The MPS is still missing the module’s unique calibration coefficients. A copy of these coefficients is supplied on a CD with the module. Using the same procedure as Step #6, drag-and-drop the module’s calibration coefficients (CAL.CFG file). Cycling power again will pull these calibration coefficients into the working RAM. If the calibration coefficients cannot be located, please contact Scanivalve Corp and provide your module’s exact model and serial number.

[INTENTIONALLY LEFT BLANK]

[INTENTIONALLY LEFT BLANK]

APPENDIX

APPENDIX A: INTRA-FRAME A/D TIMING

The rate setting controls the rate that each frame is sampled. The output rate may vary slightly based on the network loading and latencies. This section describes the timing of each pressure channel within a frame.

The hardware consists of 16 A/D chips, with each chip having 4 inputs. This arrangement produces 64 channels of pressure data. The 4 A/D channels are read sequentially. The relationship between the A/D chips, A/D channels and the pressure channels are shown in the map below:

A/D	CHIP 1	CHIP 2	CHIP 3	CHIP 4	CHIP 5	CHIP 6	CHIP 7	CHIP 8
CHAN 1	1	5	9	13	17	21	25	29
CHAN 2	2	6	10	14	18	22	26	30
CHAN 3	3	7	11	15	19	23	27	31
CHAN 4	4	8	12	16	20	24	28	32

A/D	CHIP 9	CHIP 10	CHIP 11	CHIP 12	CHIP 13	CHIP 14	CHIP 15	CHIP 16
CHAN 1	36	40	44	48	52	56	60	64
CHAN 2	35	39	43	47	51	55	59	63
CHAN 3	34	38	42	46	50	54	58	62
CHAN 4	33	37	41	45	49	53	57	61

The following table shows the timing of the start of each pressure channel A/D sample relative to the start of each frame:

TIME	CHAN	TIME (NOTE 1)	CHAN	TIME (NOTE 1 & 2)	CHAN	TIME (NOTE 1)	CHAN
0.0us	1	P/4 + 0us	2	2P/4 + 0us	3	3P/4 + 0us	4
4.4us	5	P/4 + 4.4us	6	2P/4 + 4.4us	7	3P/4 + 4.4us	8
8.8us	9	P/4 + 8.8us	10	2P/4 + 8.8us	11	3P/4 + 8.8us	12
13.2us	13	P/4 + 13.2us	14	2P/4 + 13.2us	15	3P/4 + 13.2us	16
17.6us	17	P/4 + 17.6us	18	2P/4 + 17.6us	19	3P/4 + 17.6us	20
22.0us	21	P/4 + 22us	22	2P/4 + 22us	23	3P/4 + 22us	24
26.4us	25	P/4 + 26.4us	26	2P/4 + 26.4us	27	3P/4 + 26.4us	28
30.8us	29	P/4 + 30.8us	30	2P/4 + 30.8us	31	3P/4 + 30.8us	32
35.2us	36	P/4 + 35.2us	35	2P/4 + 35.2us	34	3P/4 + 35.2us	33
39.6us	40	P/4 + 39.6us	39	2P/4 + 39.6us	38	3P/4 + 39.6us	37
44.0us	44	P/4 + 44.0us	43	2P/4 + 44.0us	42	3P/4 + 44.0us	41
48.4us	48	P/4 + 48.4us	47	2P/4 + 48.4us	46	3P/4 + 48.4us	45
52.8us	52	P/4 + 52.8us	51	2P/4 + 52.8us	50	3P/4 + 52.8us	49
57.2us	56	P/4 + 57.2us	55	2P/4 + 57.2us	54	3P/4 + 57.2us	53
61.6us	60	P/4 + 61.6us	59	2P/4 + 61.6us	58	3P/4 + 61.6us	57
66.0us	64	P/4 + 66.0us	63	2P/4 + 66.0us	62	3P/4 + 66.0us	61

Notes:

1. P = 1 / Rate
2. 8.4us is added to all times in this column due to temperature read during scan.

APPENDIX B - INTERNAL SOCKET ALLOCATION

The following is the socket allocation for the MPS 4000:

SOCKET ALLOCATION	
SOCKET	ALLOCATION
0	PHY Control Frames
1	FTP Server Control (Primary) – when not scanning FTP Client – when scanning Calibration Client MPS Master Telnet Client to Slave
2	Data Server (Primary) FTP Data Client FTP Data Server MPS Master Oven Client DHCP Communication
3	PTP Event 319
4	PTP General 320
5	Web Server
6 *	Multicast Server (Default, set with SVRSEL set to 1) Binary Server (Set with SVRSEL set to 2)
7	Telnet Server

* Note: Socket 6 action is selected with the SVRSEL variable. Changing this variable's setting requires a reboot.

APPENDIX C - LABVIEW CLIENT EXAMPLE

Below is an example LabVIEW VI to show how to interface with the MPS 4000 binary server. The frame structure to the left connects to the MPS4000 binary server with port 503 and prompts the user to enter a file name for the text file that will hold the scan data in CSV format.

The next frame structure to the right creates a queue and starts the scan by writing a 1 to the binary server. A queue is used between the producer loop and the consumer loop to handle differences between sampling rates and writing rates.


The next frame structure to the right supports the produce loop and the consumer loop.

The producer loop reads the data from the MPS4000 binary server. The MPS data is formatted such that the VI cast function can convert it directly to 66 single elements per frame or 264 bytes. This example runs the produce loop at 1ms, but this could be run at a slower period when the rate of the MPS is slower.

Multiple frames of data from the MPS may be read at once, however this will cause multiple frames of data in one array. For instance if 4 frames of data are desired with one read, the number of bytes to read is (4 X 264) or 1056 bytes. This allows the period of the producer loop to be reduced but the spreadsheet will then have 4 frames of data in one row. Additional LabVIEW functions, or spreadsheet post processing functions, are required to split out each data frame when multiple frames are read.

In tests with the produce loop set at a period of 1ms, and the read size set equal to 1 frame, data rates to 850 Hz are possible.

The consumer loop reads data from the queue, displays it to the front panel, converts the data to a string formatted for a CSV file, and writes the data to a CSV file. When an error is received from the queue, because the queue was closed from the producer loop, the consumer loop exits and closes the CSV file.

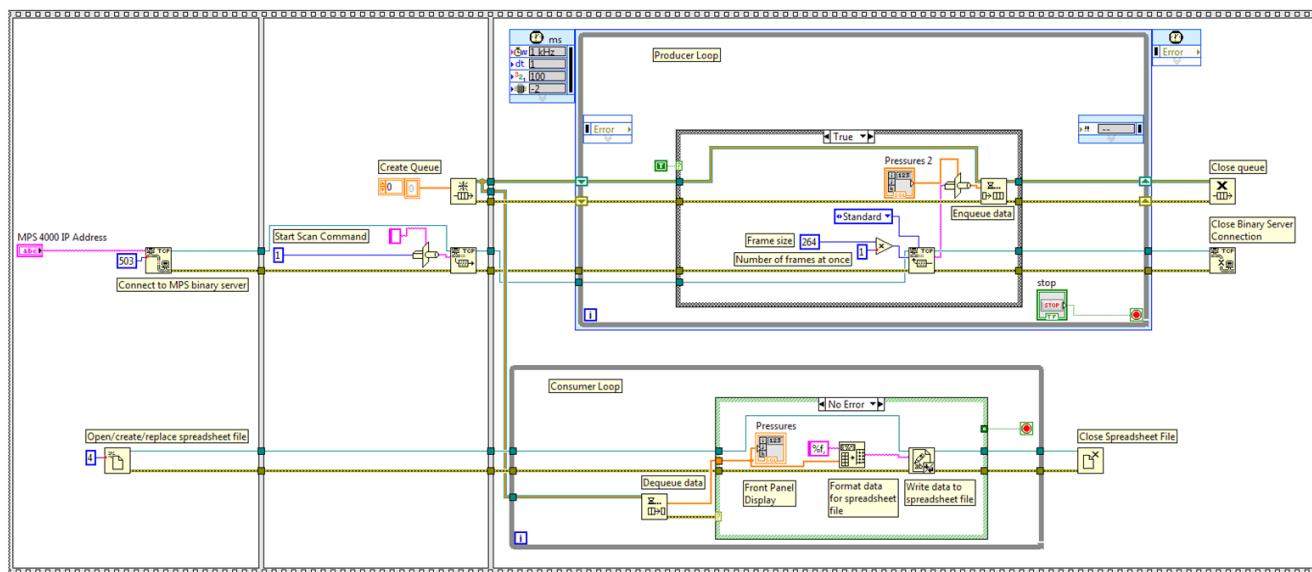

APPENDIX D: 1588 SYNCHRONIZED SCAN START

The 1588 protocol allows multiple MPS4000 units to start scanning at a preset time while keeping their data sample time synchronized. The scan start preset time, SSD and SST, could be a future time or a time past. The following tables show the time line for both cases. Please note that the frame timestamp is always relative to the start of scan time.

The following table describes a start scan with two MPS4000 units with internal trigger when the setting of SSD and SST are in the future.

SCAN START WITH FUTURE START TIME FOR 2 MPS4264 MODULES, INTERNAL TRIGGER							
SSD SET TO 6/1/2020							
SST SET TO 12:00:00							
RATE SET TO 1Hz							
1588 DATE AND TIME	EVENT	MPS A FRAME	MPS B FRAME	MPS A FRAME TIME STAMP	MPS B FRAME TIME STAMP	MPS A EXT TRIG TIME STAMP	MPS B EXT TRIG TIME STAMP
6/1/2020 11:59:58.00							
6/1/2020 11:59:58.50	Scan to MPS A						
6/1/2020 11:59:59.00							
6/1/2020 11:59:59.75	Scan to MPS B						
6/1/2020 12:00:00.00							
6/1/2020 12:00:01.00		1	1	1.0	1.0	0	0
6/1/2020 12:00:02.00		2	2	2.0	2.0	0	0

The following table describes a start scan with two MPS4000 units with internal trigger when the setting of SSD and SST are in the past. The MPS always starts on a period (1/rate) interval from the setting of the SSD and SST.

SCAN START WITH PAST START TIME FOR 2 MPS4264 MODULES, INTERNAL TRIGGER								
SSD SET TO 6/1/2020								
SST SET TO 11:59:58								
RATE SET TO 1Hz								
1588 DATE AND TIME	EVENT	MPS A FRAME	MPS B FRAME	MPS A FRAME TIME STAMP	MPS B FRAME TIME STAMP	MPS A EXT TRIG TIME STAMP	MPS B EXT TRIG TIME STAMP	MPS B EXT TRIG TIME STAMP
6/1/2020 11:59:58.00								
6/1/2020 11:59:59.00								
6/1/2020 11:59:59.75	Scan to MPS A							
6/1/2020 12:00:00.00								
6/1/2020 12:00:00.25	Scan to MPS B							
6/1/2020 12:00:01.00		1		1.0		0		
6/1/2020 12:00:02.00		2	1	2.0	1.0	0	0	
6/1/2020 12:00:03.00		3	2	3.0	2.0	0	0	

The following table describes a start scan with two MPS4000 units with external trigger when the setting of SSD and SST are in the past. The MPS always starts on a period (1/rate) interval from the setting of the SSD and SST. The external trigger time stamp indicates when the external trigger occurred in 1588 time.


SCAN START WITH PAST START TIME FOR 2 MPS4264 MODULES, EXTERNAL TRIGGER								
SSD SET TO 6/1/2020								
SST SET TO 11:59:58								
RATE SET TO 1Hz								
1588 DATE AND TIME	EVENT	MPS A FRAME	MPS B FRAME	MPS A FRAME TIME STAMP	MPS B FRAME TIME STAMP	MPS A EXT TRIG TIME STAMP	MPS B EXT TRIG TIME STAMP	MPS B EXT TRIG TIME STAMP
6/1/2020 11:59:58.00								
6/1/2020 11:59:59.00								
6/1/2020 11:59:59.75	Scan to MPS A							
6/1/2020 12:00:00.00								
6/1/2020 12:00:00.25	Scan to MPS B							
6/1/2020 12:00:01.00								
6/1/2020 12:00:02.00								
6/1/2020 12:00:02.57	External Trig							
6/1/2020 12:00:03.00		1	1	2.0	2.0	6/1/2020 12:00:02.57	6/1/2020 12:00:02.57	
6/1/2020 12:00:04.00								
6/1/2020 12:00:04.29	External Trig							
6/1/2020 12:00:05.00		2	2	4.0	4.0	6/1/2020 12:00:04.29	6/1/2020 12:00:04.29	

APPENDIX E - ENGINEERING UNIT CONVERSION CONSTANTS


UNITSCAN Setting	Engineering Unit	PSI to EU 1 psi =	EU to PSI 1 EU =
ATM	Atmospheres	0.068046 A	14.6960 psi
BAR	Bars	0.068947 b	14.5039 psi
CMHG	Centimeter of Mercury	5.17149 cmHg	0.193368 psi
CMH2O	Centimeter of Water	70.308 cmH2O	0.014223 psi
DECIBAR	Decibar	0.68947 db	1.4504 psi
FTH2O	Foot of Water	2.3067 ftH2O	0.43352 psi
GCM2	Gram per square Centimeter	70.306 g/cm2	0.014224 psi
INHG	Inch of Mercury @ 0C	2.0360 inHg	0.491159 psi
INH2O	Inch of Water @ 4C	27.680 inH2O	0.036127 psi
KGCM2	Kilogram per square Centimeter	0.0703070 kg/cm2	14.2235 psi
KGM2	Kilogram per square Meter	703.069 kg/m2	0.0014223 psi
KIPIN2	kips per square inch(ksi)	0.001 kip/in2	1000.0 psi
KNM2	Kilonewton per square Meter	6.89476 kN/m2	0.145038 psi
KPA	Kilopascal	6.89476 kPa	0.145038 psi
MBAR	Millibar	68.947 mb	0.014504 psi
MH2O	Meter of Water	0.70309 mH2O	1.42229 psi
MMHG	Millimeter of Mercury	51.7149 mmHg	0.0193368 psi
MPA	Megapascal	0.00689476 Mpa	145.038 psi
NCM2	Newton per square Centimeter	0.689476 N/cm2	1.45038 psi
NM2	Newton per square Meter	6894.759766 N/m2	0.000145038 psi
OZFT2	Ounce per square Foot	2304.00 oz/ft2	0.000434028 psi
OZIN2	Ounce per square Inch	16.00 in/ft2	0.062500 psi
PA	Pascal	6894.759766 Pa	0.000145038 psi
PSF	Pound per square Foot	144.00 lb/ft2	0.00694444 psi
TORR	Torr	51.714901 T	0.0193368 psi



APPENDIX F - CABLE DRAWINGS AND PINOUTS



MPS PREMIUM POWER/SERIAL/trigger CABLE

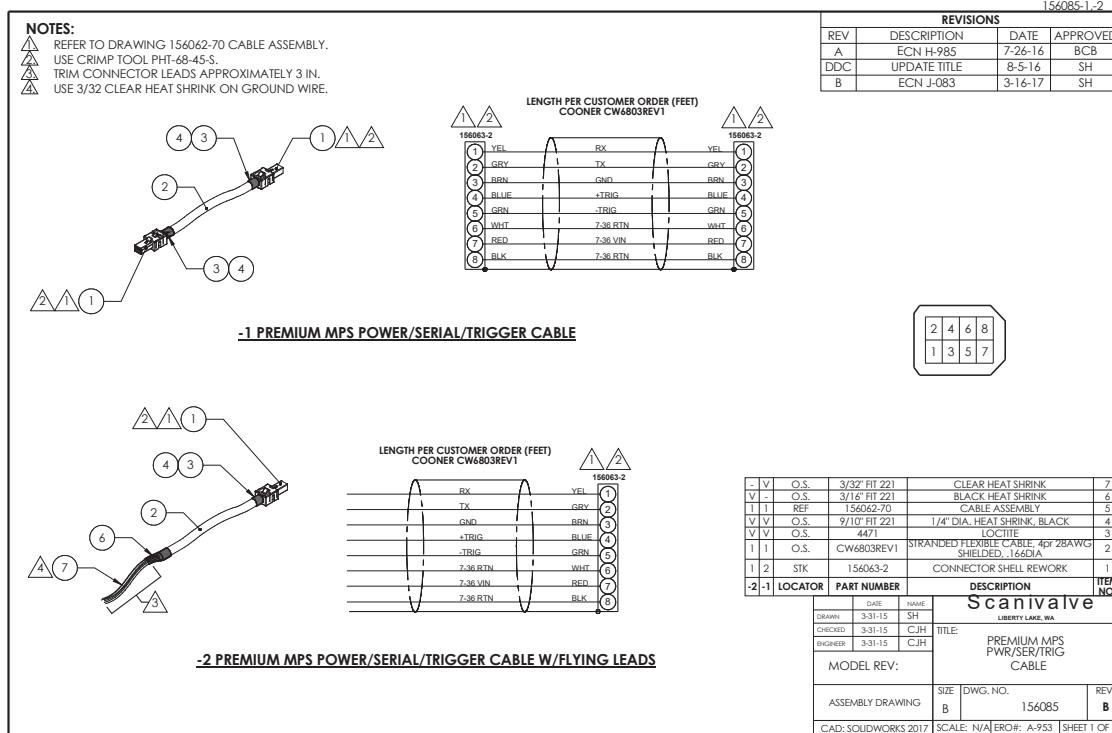

MPS POWER/trigger CABLE WITH FLYING LEADS

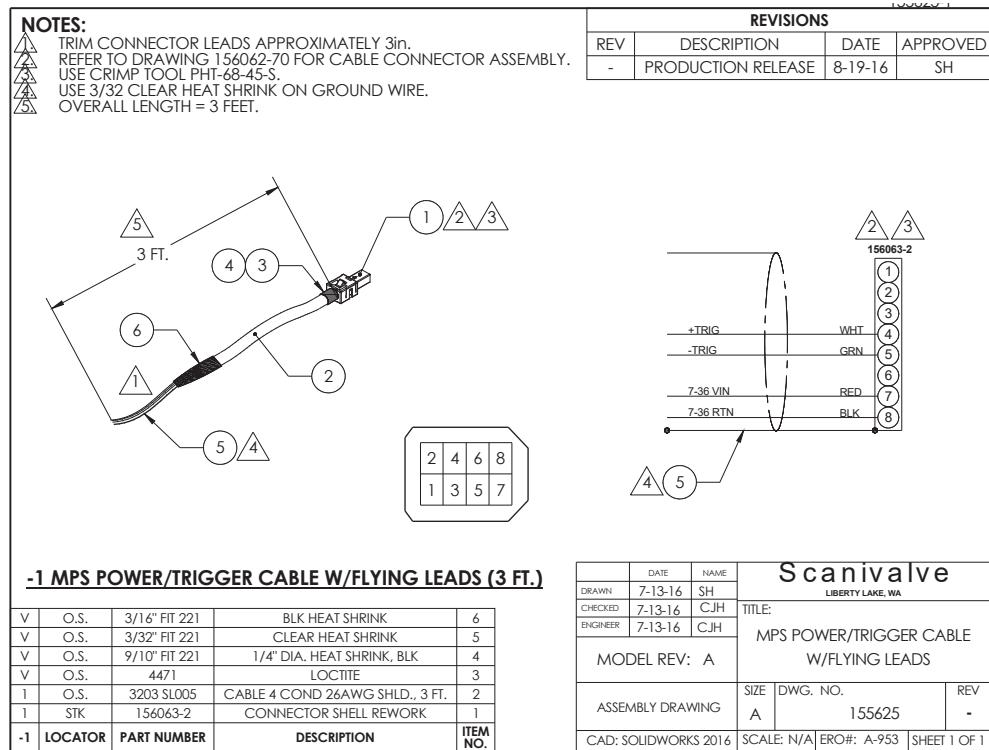

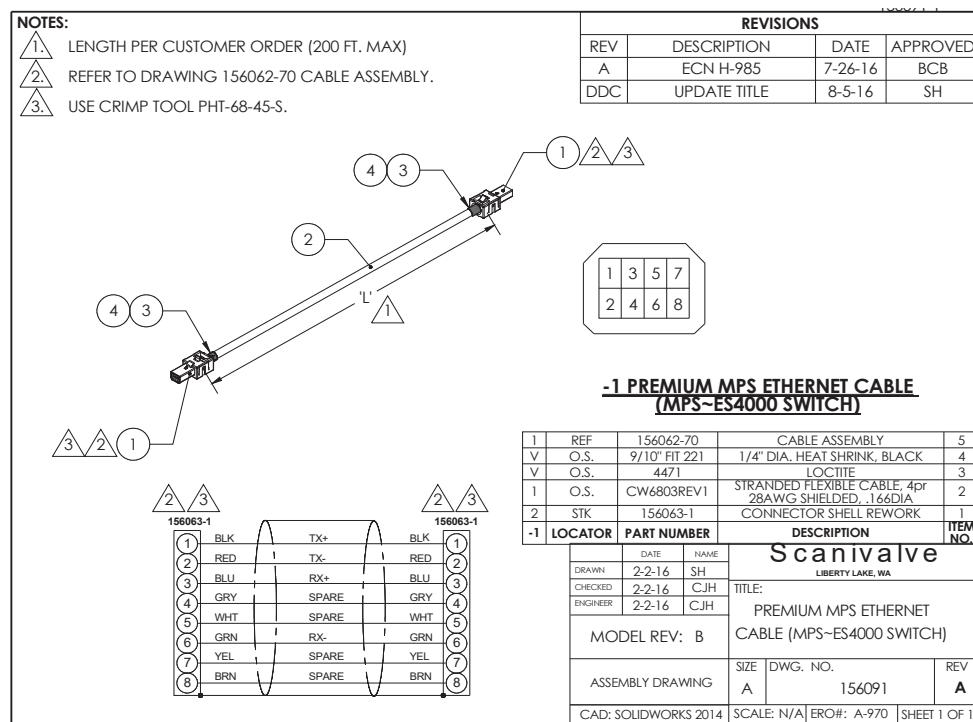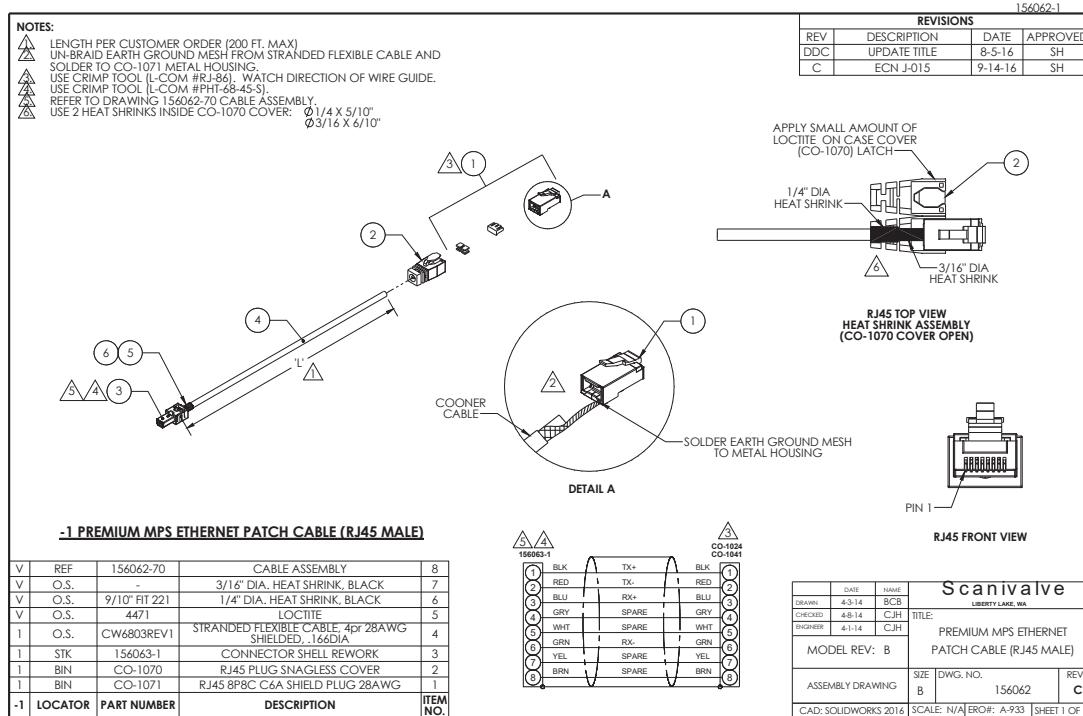MPS PREMIUM ETHERNET CABLE**MPS TO RJ45 ETHERNET CABLE**

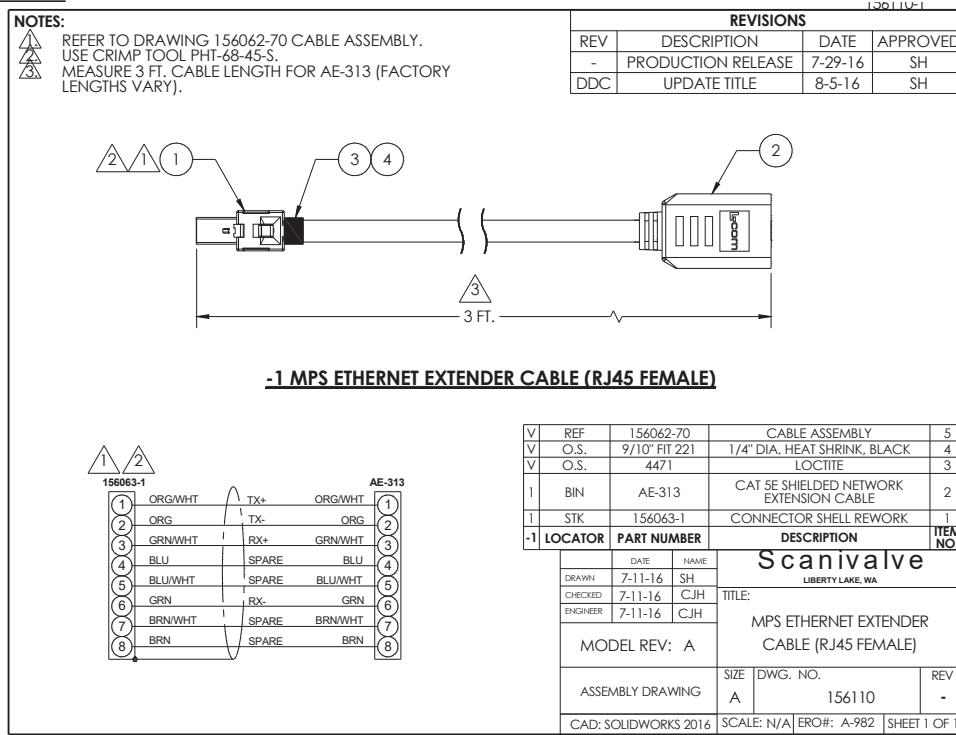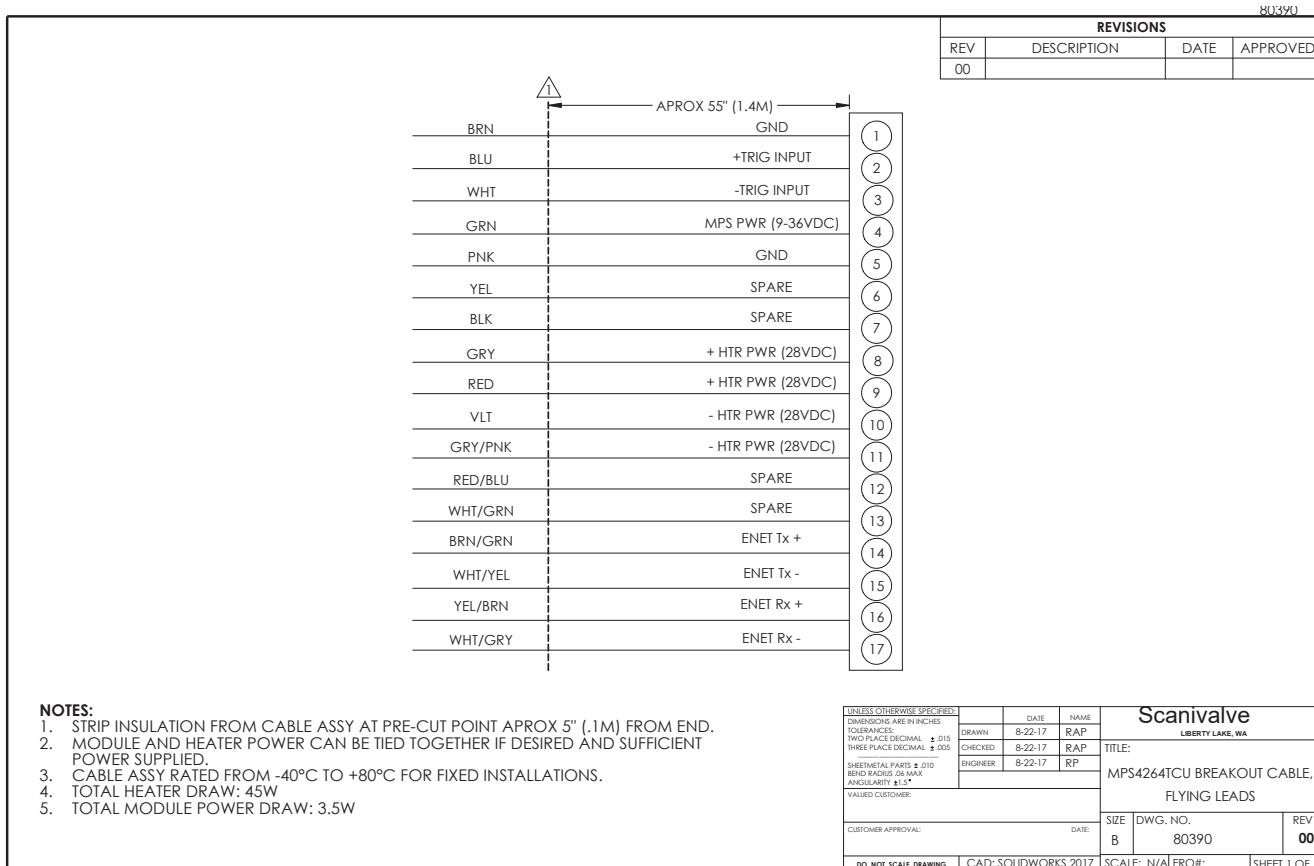MPS TO RJ45 FEMALE**MPS4264TCU CABLE**

APPENDIX G - SOFTWARE REVISION LOG

Version 1.00 - Never released

Version 2.00 - Initial release (8/24/2016)

Version 2.01 - Improved temperature mapping of sensors. Improved reliability of calibration via SPC4000. (2/9/17)

Version 2.02 - Added median filtering to calibration. Added auto multi-unit validation for users with an SPC4000.(2/24/17)

Version 2.03 - Resolved a bug that didn't allow the current temperature plane to be read before a CALZ. (3/6/17)

Version 2.04 - Added automatic and manual field validation capabilities for single and multiple units. Allowed a Telnet connection to be established while in the webserver. (4/28/2017)

Version 2.05 - Added support for multiple calibration ovens. Added the ability for an MPS to act as a PTP Master. Resolved an issue with the scan rate and time with the binary packet over the binary server. Resolved a bug that altered the temperature averaging of the Binary LabVIEW packet. Resolved a bug that caused errors when scanning multiple units via FTP. (11/9/2017)

Version 2.06 - Removed SAVE C and SAVE T from the SAVE command. Resolved an issue with SET FORMAT. Added new error routing. Added DINT and EINT around clock update. Added integral to servo control loop. Resolved PTP resest issue. (9/5/2018)

Version 2.07 - Improved reliability of FTP data collection at high scan rates. (1/11/19)

[INTENTIONALLY LEFT BLANK]

Scanivalve

**1722 N. MADSON ST.
LIBERTY LAKE, WA 99019**

**PHONE: 1-800-935-5151
1-509-891-9970
FAX: 1-509-891-9481**

**SCANCO@SCANIVALVE.COM
WWW.SCANIVALVE.COM**

**MPS4264 HARDWARE & SOFTWARE MANUAL
JANUARY 11TH, 2019**